

**LABORATORY
NOTEBOOK**

Nº 1044

medpointe

**PLAINTIFFS'
TRIAL EXHIBIT**
PTX0142

exhibitsticker.com

CIP2061
Argentum Pharmaceuticals LLC v. Cipla Ltd.
IPR2017-00807

HIGHLY CONFIDENTIAL
-SUBJECT TO STIPULATED PROTECTIVE ORDER

MEDA_APTX03470948

NOTEBOOK NO. R 1044
ISSUED TO John D'Acosta
ON August 5 2005
DEPARTMENT Product & Process Development
RETURNED November 29 2007

LABORATORY NOTEBOOK ERROR CODES

CE CALCULATION ERROR
ER ENTRY ERROR
GE GRAMMATICAL ERROR
LE LATE ENTRY
RE ROUNDING ERROR
SE SPELLING ERROR
TE TRANSCRIPTION ERROR
WE WORDING ERROR

SCIENTIFIC NOTEBOOK COMPANY
2831 LAWRENCE AVENUE
STEVENSVILLE, MICHIGAN 49127
(800) 537-3028 - <http://www.snco.com>

INSTRUCTIONS

1. The primary purpose of this notebook is to protect your and the Company's Patent-Rights by keeping records of all original work in a form acceptable as evidence if any legal conflict arises.

2. When starting a page, enter the title, project number, and book number.

• Use ink for permanence -- avoid pencil.

• Record your work as you progress, including any spur-of-the-moment ideas which may be developed later.

• Avoid making notes on loose paper to be recopied.

• Record your work in such a manner that a co-worker can continue from where you stop. You might be ill and to protect your priority it could be urgent that the work continue while you are absent.

3. Give a complete account of your experiments and the results, both positive and negative, including your observations. Record all diagrams, layouts, plans, procedures, new ideas, or anything pertinent to your work including the details of any discussions with suppliers, or other people outside the Company.

• Do not try to erase any incorrect entries; draw lines deleting them, note the corrections, sign and date the changes. This extra care is worthwhile because of the necessity of original data to prove priority of new discoveries.

4. After entering your data, sign and date the entries.

• Explain your work to at least two wit-

nesses who are not co-inventors, and have them sign and date the pages in the place provided.

• Record the names of operators and witnesses present during any demonstration and have at least two witnesses sign the page. If no witnesses are present during an experiment of importance, repeat it in the presence of two witnesses.

5. Since computer programs can be patented these instructions apply to the development of computer software. In this case a description of the structure and operation of the program should be recorded in the notebook, together with a basic flow diagram which illustrates the essential features of the program. In the course of developing the code, the number of lines of code written each day should be recorded in the notebook, together with a statement of the portion of the flow diagram to which the section of code is directed.

6. This notebook and its contents are the exclusive property of the Company. It is confidential and the contents are not to be disclosed to anyone unless authorized by the Company. You must return it when completed, upon request, or upon termination of employment. It should be kept in a protected place. **If loss occurs, notify your supervisor immediately, and make a written report describing the circumstances of the loss.**

Copyright Harold Gallop

Form 590

Table of Contents

	Page	DATE
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-146A WET GRANULATION	1	08/12/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-146B WET GRANULATION	2	08/12/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-148A WET GRANULATION	3	08/17/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-148B WET GRANULATION	4	08/17/05
SOMA COMBO: EXCIPIENT COMPATIBILITY STUDY 2 WEEK PULL @ 50°C	5	08/19/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-150A WET GRANULATION	6	08/19/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-150B WET GRANULATION	7	08/19/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-152A WET GRANULATION	8	08/25/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-152B WET GRANULATION	9	08/25/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-154 WET GRANULATION	10	08/29/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-156 WET GRANULATION	11	08/29/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-160 WET GRANULATION	12	08/31/05
TUSSI-12 DS SUSPENSION: SAMPLES RECEIVED FROM DECATUR. FOUR-ABILITY OBSERVATIONS & VISCOSITY TESTING	13	09/01/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-162 WET GRANULATION	14	08/31/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-164 WET GRANULATION	15	09/06/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-168 WET GRANULATION	16	09/07/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-170 WET GRANULATION	17	09/08/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-172 WET GRANULATION	18	09/09/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-174 WET GRANULATION	19	09/19/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-176 WET GRANULATION	20	09/20/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-178 WET GRANULATION	21	09/20/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-180 WET GRANULATION	22	09/21/05
VISCOSITY EXPERIMENT: PECTIN LOT #0000001174	23	09/21/05

SCIENTIFIC NOTEBOOK CO.

PTX0142-00003

MEDA_APTX03470950

Table of Contents

Page

DATE

STABILITY: SOMA SR #1009-159, 160, 162, 164, 168, 170, 172, 176, & 178 } @ 20kN PLACED ON 25/10 & 40/75 STABILITY.	24	09/22/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-182 WET GRANULATION	25	09/22/05
PHYSICAL TESTING: SOMA 700mg SR TABLETS #1009-184 WET GRANULATION	26	09/23/05
STABILITY TESTING: TUSSI AM SUSPENSION #1027-21 @ 1 MONTH	27	09/23/05
TUSSI AM REMAKE OF #1030-3 TO COMPARE VISCOSITY RESULTS (PECTIN LOT #000001174)	28	10/05/05
TUSSI AM REMAKE OF #1030-3 AND #1044-28 TO COMPARE VISCOSITY RESULTS (PECTIN LOT #000002012)	30	10/26/05
PHYSICAL TESTING: SOMA SR #1048-2, 4, 6, 8 WET GRANULATIONS	32	10/06/05
TAP DENSITY TESTING: SOMA SR #1048-2, 4, 6, 8 WET GRANULATIONS	33	10/09/05
PHYSICAL TESTING: SOMA SR #1048-10, 12, 14, 16 WET GRANULATIONS	34	10/09/05
TAP DENSITY TESTING: SOMA SR #1048-10, 12, 14, 16 WET GRANULATIONS	35	10/10/05
PARTICLE SIZE DISTRIBUTION: SOMA SR #1048-2, 4, 6, 8	36	10/11/05
SOLUBILITY EXPERIMENT: METHOCEL K100-LV-OR	38	10/11/05
PHYSICAL TESTING: SOMA SR #1048-18 & 20 WET GRANULATIONS	39	10/12/05
PARTICLE SIZE DISTRIBUTION: SOMA SR #1048-10, 12, 14, 16	40	10/12/05
PHYSICAL TESTING: SOMA SR #1048-22 WET GRANULATION	42	10/13/05
PHYSICAL TESTING: SOMA SR #1048-24 WET GRANULATION	43	10/17/05
PHYSICAL TESTING: SOMA SR #1048-26 & 28 WET GRANULATIONS	44	10/19/05
PHYSICAL TESTING: SOMA SR #1048-28 & 30 WET GRANULATIONS	45	10/20/05
PHYSICAL TESTING: SOMA SR #1048-32 WET GRANULATIONS	46	10/21/05
SOLUBILITY EXPERIMENT: METHOCEL E5-LV	47	10/21/05
WEIGHT & HARDNESS VARIATION: SOMA SR #1048-36	48	10/25/05
PHYSICAL TESTING: SOMA SR #1048-36 WET GRANULATION (SCALE-UP 3kg)	50	10/26/05
TAP DENSITY & PARTICLE SIZE DISTRIBUTION: SOMA SR #1048-36	50	10/26/05

SCIENTIFIC NOTEBOOK CO.

Table of Contents

Page

DATE

STABILITY TESTING: SOMA 700mg SR #1048-24, 6, 8 @ 4 WEEKS	51	11/02/05
TUSSI AM REMAKE OF #1048-28 WITH SAME LOT OF PECTIN TO ATTEMPT TO SHOW VISCOSITY REPEATABILITY.	52	11/02/05
TUSSI AM SUSPENSION DECREASED HYDRATION TIME IN VEEGUM STEP.	54	11/04/05
TUSSI AM SUSPENSION INCREASED HYDRATION TIME IN VEEGUM STEP.	56	11/04/05
STABILITY TESTING: TUSSI AM #1044-28 @ 1 MONTH. WAS NOT PERFORMED. BATCH MADE FOR INIT. REQS.	58	11/05/05
PHYSICAL TESTING: SOMA SR #1048-39, 40, 41 WET GRANULATION	59	11/10/05
TUSSI AM SUSPENSION DECREASED HYDRATION TIME IN VEEGUM STEP.	60	11/08/05
TUSSI AM SUSPENSION INCREASED HYDRATION TIME IN VEEGUM STEP.	62	11/09/05
STABILITY TESTING: TUSSI AM #1027-21 @ 3 MONTHS. WAS NOT PERFORMED DUE TO CHAMBER DOWNTIME :: FULL DATE	64	11/09/05
ASTELIN NASAL SPRAY FORMULATION #03-33 REMAKE WITH TALIN @ REPLACING SUCRALOSE @ 0.150% SOMA SR #1048-10, 12, 14, 16	65	11/17/05
STABILITY TESTING: @ 4 WEEKS	66	11/14/05
PHYSICAL TESTING: SOMA SR #1048-42, 44, 46 WET GRANULATIONS	68	11/16/05
PHYSICAL TESTING: SOMA SR #1038-15 DECATOR WET GRANULATION @ 10, 15, 20, 25, 30 SCU + BULK	69	11/16/05
PARTICLE SIZE DISTRIBUTION & TAP DENSITY: SOMA SR #1038-15 GRAN.	70	11/16/05
PHYSICAL TESTING: SOMA SR #1048-48 & 50 WET GRANULATIONS	71	11/22/05
PHYSICAL TESTING: SOMA SR #1048-53, 54, 55 WET GRANULATIONS	72	11/22/05
PHYSICAL TESTING: SOMA SR #1048-56, 58, 60, 62 WET GRANULATIONS	73	12/01/05
PHYSICAL TESTING: SOMA SR #1048-66, 68, 70, 72, 74 WET GRANULATIONS	74	01/05/06
STABILITY: SOMA SR #1038-16, 17, 18 @ 20kN/SCU PLACED @ 25/10 44/75, 25/10	75	01/05/06
PHYSICAL TESTING: SOMA SR #1048-76, 78, 80, 82, 84 WET GRANULATIONS	76	01/05/06
PHYSICAL TESTING: SOMA SR #1048-86, 88, 90, 92, 94 WET GRANULATIONS	77	01/10/06
PHYSICAL TESTING: SOMA SR #1048-96, 98, 100, 102, 104, 108 WET GRANULATIONS	78	01/11/06
PHYSICAL TESTING: SOMA SR #1038-16 WET GRANULATION	79	01/11/06

SCIENTIFIC NOTEBOOK CO.

Table of Contents

Page

DATE

PHYSICAL TESTING: SOMA SR # 1038-17 WET GRANULATION	80	01/12/06
PHYSICAL TESTING: SOMA SR # 1038-18 WET GRANULATION	81	01/12/06
NY-TANNIC TABLETS RECEIVED FROM DR. SACKS, 10 TABLET BULK WT. 25 TABLETS SENT TO DECATUR	82	01/24/06
STABILITY: SOMA SR 700mg 1038-15-18 OPEN DISH / 1038-19-22 W/ W/ DESICCANT / 1048-10, 12 WITHOUT DESICCANT ONLY / 1048-82, 88, 108-120 W/ W/	83	01/30/06
STABILITY: SOMA SR 700mg STABILITY CHAMBERS BACK ON.	84	02/07/06
STABILITY: SOMA SR 700mg BATCH CHANGES AS PER BEN JOHNS.	85	02/08/06
PHYSICAL TESTING: SOMA SR # 1048-122 WET GRANULATION	86	02/09/06
PHYSICAL TESTING: SOMA SR # 1048-124 WET GRANULATION	87	02/10/06
PHYSICAL TESTING: SOMA SR # 1048-126 WET GRANULATION	88	02/13/06
STABILITY: SOMA SR 700mg # 1038-15, 16, 17, 18 ONE WEEK OPEN DISH	89	02/14/06
STABILITY: SOMA SR 700mg # 1038-18 RETURNED TO OPEN DISH STUDY FOR @ WK. FULL.	90	02/17/06
PHYSICAL TESTING: SOMA SR # 1048-112, 114, 116, 118, 120	91	02/20/06
PHYSICAL TESTING: SOMA SR # 1038-23	92	02/20/06
NY-TANNIC TABLETS RETURNING COMMERCIAL BOTTLE, RETAINING 50 TABLETS & KEEPING THEM IN RH. 134	94	05/20/06
EXPERIMENT: ASTELIN + STEROID COMPAT. Y STUDY BETWEEN FLUTICASONES & VP3	95	03/20/06
PHYSICAL TESTING: SOMA SR # 1048-12 & 82 COMPARE INITIAL WITH # 1038-12	96	03/21/06
PHYSICAL TESTING: SOMA SR # 1038-19, 20, 21, 22 INITIAL PERFORMED BY A-F	97	03/21/06
PHYSICAL TESTING: SOMA SR # 1038-24, 25, 26 INITIAL PERFORMED BY A-F	98	03/21/06
ASTELIN 03-33 REGISTRATION BATCHES 03-33-03c, 03-33-05s, 03-33-04s, 03-33p-03s BATCH NOTES	99	03/21/06
TUSSI-ORGANIDIN NR LIQUID CODEINE PLACEBO OF TUSSI-ORGANIDIN NR LIQUID FOR H-MORTKO FOR ANALYTICAL TESTING	102	04/05/06
TUSSI-ORGANIDIN NR LIQUID CODEINE PLACEBO OF T.O. NR LIQ. FOR H-MORTKO FOR ANALYTICAL TESTING - REDO	103	04/06/06
ASTELIN 03-33 REGISTRATION BATCHES FILTRATION STUDY RE-BONE AND GRAPHED.	104	03/24/06
ASTELIN # 1027-17 FREEZE/THAW @ CYCLES.	107	03/27/06

SCIENTIFIC NOTEBOOK CO.

Table of Contents

Page

DATE

ASTELIN REGISTRATION BATCHES KALIDAS KALE'S POST-MANUFACTURING OBSERVATION REPORT	108	03/27/06
ASTELIN 03-33 2.0 LITER BATCH FOR C-STEWART FOR LEECHABLE/EXTRACTABLE TESTING	110	04/07/06
ASTELIN 03-33p 2.0 LITER BATCH (PLACEBO) FOR L/E TESTING.	111	04/07/06
OSMOLARITY: ADJUSTMENT OF OPTIVAR WITH SORBITOL SOLUTION OPTIVAR FEASIBILITY # 4, LOT # 02306A	112	04/14/06
ASTELIN + STEROID EXPERIMENT ADDITION OF 0.1% AZELASTINE HCL TO FLUTICASONES N.S. 50mg	114	04/24/06
STABILITY TESTING: SOMA SR # 1048-10, 12, 14, 66, & 82 @ 2 MONTHS IN 70°C / 75% RH	115	04/25/06
ASTELIN + STEROID EXPERIMENT FLUTICASONES SUPERNATANT MICROSCOPY. ADDITIONAL OF AZELASTINE HCL, MICRO-PICS	116	04/25/06
ASTELIN + STEROID EXPERIMENT FLUTICASONES SUPERNATANT, ADDITION OF AZELASTINE HCL, OVERNIGHT MIX.	117	04/26/06
ASTELIN + STEROID EXPERIMENT HEATING AND CONTINUED MIXING OF # 116 SAMPLE. LEFT OVERNIGHT TO SETTLE.	118	04/27/06
ASTELIN + STEROID EXPERIMENT PREPARATION OF FLUTICASONES SUPERNATANT SAMPLE FOR ANALYTICAL + FILTERING OF SAMPLE p. 117.	119	04/28/06
PHYSICAL TESTING: SOMA 700mg SR TABLETS 27-18 & 19-01c CLINICAL BATCHES	120	04/28/06
SOLUBILITY EXPERIMENT: AZELASTINE HCL IN PEG 400	121	05/01/06
OSMOLARITY TESTING: OPTIVAR/SORBITOL SOLUTION SAMPLES PREPARED @ 5%, 7.0%, 7.5% TO TEST FOR SOLUBILITY DIFFERENCES	122	05/02/06
OSMOLARITY TESTING: OPTIVAR/SORBITOL SOLUTION PREPARATION OF (3) OPTIVAR BATCHES, WITH VARYING SORBITOL SOL. (70%)	123	05/03/06
OPTIVAR: CONCENTRATIONS, FOR OSMOLARITY TESTING.	124	05/03/06
ASTELIN NASAL SPRAY # 03-33 FOR A.E.T. TESTING B.A.C. @ 80%	126	05/11/06
ASTELIN NASAL SPRAY # 03-33 FOR A.E.T. TESTING B.A.C. @ 50%	127	05/12/06
ASTELIN NASAL SPRAY # 03-33 FOR A.E.T. TESTING B.A.C. @ 25%	128	05/15/06
ASTELIN NASAL SPRAY # D3-33pp FOR A.E.T. TESTING PRESERV. PLACEBO, B.A.C. = 0	129	05/16/06
ASTELIN NASAL SPRAY # 03-36 FOR A.E.T. TESTING B.A.C. @ 80%	130	05/11/06
ASTELIN NASAL SPRAY # 03-36 FOR A.E.T. TESTING B.A.C. @ 50%	131	05/12/06
ASTELIN NASAL SPRAY # 03-36 FOR A.E.T. TESTING B.A.C. @ 25%	132	05/15/06
ASTELIN NASAL SPRAY # 03-36pp FOR A.E.T. TESTING B.A.C. = 0, PRESERV. PLACEBO	133	05/16/06

SCIENTIFIC NOTEBOOK CO.

Table of Contents	Page
OSMOLARITY TESTING: NEW STANDARDS ARRIVE: WILL ATTEMPT TO CALIBRATE AND RE-TEST #1044-124 a, b, c.	134 05/21/06
OPTIVAR UNIT DOSE: DELIVERY OF (28) ^{Lot # 01406A} CASES FROM CARDINAL HEALTH	136 06/07/06
OPTIVAR UNIT DOSE OSMOLARITY TESTING OF LOT # 01406A	137 06/07/06
AZELASTINE HCL + AVICEL RC591 EXPERIMENTS	138 06/08/06
ASTELIN NASAL SPRAY #03-33 ^{REMAKE OF #1044-129 B/C} _{ANALYTICAL FOUND BAC IN PRODUCT}	140 06/13/06
OSMOLARITY TESTING: OPTIVAR 6mL FILL, BATCH #056614 ^{FROM PATHEON}	141 06/14/06
RAW MATERIALS REC'D: FLUTICASON PROPIONATE FROM ^{HOVIONE/} _{GYMA}	142 06/15/06
OPTIVAR EXPERIMENT: END OF USP SPEC. & IMPACT ON OSMOLARITY. ^{2 BATCHES TO MIMIC SORBITOL BEING AT THE LOW}	143 06/20/06
OPTIVAR: METHYLENE BLUE DYE TEST ^{PACKAGE INTEGRITY TEST FOR OPTIVAR} _{BLOW-FILL-SEAL SAMPLES}	146 07/19/06
OPTIVAR: METHYLENE BLUE DYE ^{WE @ 10/20/06} _{TESTING OF BLUE DYE ABSORPTION w/ SPECTROPHOTOMETER IN ANALYTICAL.}	149 08/17/06
OPTIVAR: METHYLENE BLUE DYE ^{TESTING OF BLUE DYE THROUGH 0.2, 0.3, 10.4, 0.5 MICRON MICROPIPETTE}	151 08/21/06
ASTELIN + STEROID EXPERIMENTS ^{AZELASTINE & NaCMC ORDER OF ADDITION} _{EXPERIMENTS}	153 11/20/06
ASTELIN + STEROID EXPERIMENTS ^{AZELASTINE + AVICEL RC-591} _{ORDER OF ADDITION EXP.}	154 11/21/06
ASTELIN + STEROID EXPERIMENTS ^{AZELASTINE + AVICEL RC-591} _{ORDER OF ADDITION EXP.}	155 11/22/06
METHOCEL SOLUTION EXPERIMENT ^{VISCOSITY OF METHOCEL ES AND K100LVC} _{AND K100MCR SOLUTIONS}	156 02/12/07
ASTELIN 03-33 PUMP CLEANING EVAL. ^{REPORT PPD-07-07 P} _{USING 03-33-03 c}	157 03/01/07
ASTELIN 03-33 HOLD TIMES	160 03/21/07
OPTIVAR: METHYLENE BLUE DYE TEST ^{REPEAT TESTING WITH} _{0.2 & 0.5 MICRON SIZE ONLY}	161 08/27/07
N/A	

SCIENTIFIC NOTEBOOK CO.

Table of Contents	Page
N/A	

SCIENTIFIC NOTEBOOK CO.



RECEIVED NY-TAMNIC™ TABLETS FROM B. JOHNS, SAME 100 COUNT BOTTLE AS WAS LOGGED IN ON 01/24/06, AND WAS ASKED TO RETURN 50 TABLETS IN AN HDPE BOTTLE AND RETURN THE BOTTLE. THERE WERE A TOTAL OF 76 TABLETS IN THE BOTTLE AT THIS TIME. B. JOHNS RECEIVED THE NY-TAMNIC™ COMMERCIAL BOTTLE BACK FROM ME WITH ONLY 26 TABLETS. THE 50 TABLET SAMPLE WAS PLACED WITH OUR OTHER TABLET REMAINS IN RA. 1374 AND LABELED ACCORDINGLY.

26 tablets given to Janet for sending outside as per Alex. *Benjamin 3/20/06*

~~NA~~

Recorded by: *[Signature]* Date: 03/20/06
 Verified by: *[Signature]* Date: 3/20/06
 To Page No. 414

Using a commercial 17ml bottle and Vial's VP3 pump, WE ATTEMPTED TO PASS COMMERCIAL FLUTICASONIDE PROPIONATE NASAL SPRAY (LOT # C187889) THROUGH THE ASTELIN PUMPING. 7 ml = 10.98g (17ml bottle + VP3 pump, cap, & water lock) 6 pumps = 25.55g (pump assembly + bottle + suspension) NET = 14.57g

Four to six priming sprays were attempted by K. KALE. AFTER WHICH A UNIFORM PLUME WAS OBSERVED.

~~NA~~

Recorded by: *[Signature]* Date: 03/20/06
 Verified by: *[Signature]* Date: 10/25/07
 To Page No. 414

Project No. _____
Book No. _____

TITLE PHYSICAL TESTING: SOMA SR #1048-12 & 82
INITIAL U.S. 3+ MONTH 40/75. 03/21/06

From Page No. N/A

Physical Testing: Soma 700mg SR Tablets #1048-12 & 82

Upper Comp. (kN)	Target Wt. (mg)	Tablet #	Weight (g)	Hardness (SCU)	Friability
1048-12	1000	1	1.0324	23.1	Initial Wt (g) =
@		2	0.9993	17.5	10.0660
20kN		3	0.9809	15.8	Final Wt (g) =
(initial)		4	0.9899	17.2	10.0250
		5	0.9932	17.8	%Friable =
1000					
Average					
SCU					
Hardness					
1048-12	1000	1	1.0052	21.6	Initial Wt (g) =
@		2	1.0027	20.3	9.9338
20kN		3	0.9893	20.5	Final Wt (g) =
(40°C/75%RH)		4	1.0022	20.6	9.8774
		5	0.9949	20.2	%Friable =
1000					
Average					
SCU					
Hardness					
1048-82	1000	1	1.0102	28.3	Initial Wt (g) =
@		2	1.0018	25.7	10.1100
20kN		3	1.0121	27.6	Final Wt (g) =
(initial)		4	1.0163	29.4	10.0917
		5	1.0102	28.9	%Friable =
1000					
Average					
SCU					
Hardness					
1048-82	1000	1	1.0075	28.5	Initial Wt (g) =
@		2	1.0079	29.5	10.0430
20kN		3	1.0140	29.5	Final Wt (g) =
(40°C/75%RH)		4	1.0048	28.1	10.0236
		5	0.9994	25.4	%Friable =

John S. Deacony
0.0000 g
SOMA SR 700mg
1.0052 g
1.0027 g
0.9993 g
1.0022 g
0.9949 g
9.9338 g
9.8774 g
1.0075 g
1.0079 g
1.0140 g
1.0048 g
0.9994 g

1048-12
WITH
DESICCANT
1048-82

**Hardness readings in bold = tablets that capped during testing

NR	5	NR	5
XM	20.6 SC	XM	28.2 SC
SREL	2.8 %	SREL	6.0 %
SD	0.6 SC	SD	1.8 SC
XMIN	20.2 SC	XMIN	25.4 SC
XMAX	21.6 SC	XMAX	29.5 SC
5	20.2 SC	5	25.4 SC
4	20.6 SC	4	28.1 SC
ZERO		3	29.5 SC
3	20.5 SC	2	29.5 SC
ZERO		1	28.5 SC
2	20.3 SC		
1	21.6 SC		

1048-12
1048-82

N/A

To Page No. N/A

Recorded by: *John S. Deacony* Date: 03/21/06 Verified by: *QA* Date: 24 Oct 07

Project No. _____
Book No. _____

TITLE PHYSICAL TESTING: SOMA SR #1038-19, 20, 21, 22
INITIAL U.S. 3+ MONTH 40/75. 03/21/06

From Page No. N/A

Physical Testing: Soma 700mg SR Tablets #1038-19, 20, 21, & 22

Batch #	Target Wt. (mg)	Tablet #	Weight (g)	Thickness (mm)	Hardness (SCU)	Friability
1038-19	1000	1	0.9971	6.58	23.7	Initial Wt (g) =
		2	0.9947	6.60	26.2	9.9650
		3	1.0048	6.65	25.3	Final Wt (g) =
		4	0.9940	6.60	25.0	9.9410
		5	0.9951	6.58	23.8	
		6	0.9980	6.60	25.7	
		7	0.9989	6.58	26.7	
		8	0.9931	6.60	24.0	
		9	0.9912	6.57	25.3	
		10	0.9912	6.59	22.6	%Friable =
1038-20						
	900	1	0.9926	6.60	27.0	Initial Wt (g) =
		2	0.9977	6.51	24.5	9.9850
		3	0.9996	6.62	28.8	Final Wt (g) =
		4	0.9985	6.63	26.4	9.9500
		5	1.0015	6.62	28.8	
		6	0.9839	6.57	25.0	
		7	1.0011	6.65	27.0	
		8	1.0001	6.65	27.5	
		9	0.9901	6.60	25.9	
		10	1.0080	6.65	31.0	%Friable =
1038-21						
	900	1	1.0070	6.65	27.9	Initial Wt (g) =
		2	0.9983	6.58	24.3	9.9780
		3	1.0018	6.63	26.9	Final Wt (g) =
		4	0.9997	6.60	24.8	9.9370
		5	0.9914	6.57	24.0	
		6	0.9980	6.60	24.3	
		7	1.0032	6.64	29.7	
		8	0.9875	6.53	25.6	
		9	0.9980	6.59	27.6	
		10	0.9918	6.56	22.8	%Friable =
1038-22						
	900	1	1.0088	6.55	20.0	Initial Wt (g) =
		2	1.0109	6.52	25.4	10.0720
		3	1.0093	6.52	21.5	Final Wt (g) =
		4	0.9983	6.53	19.6	10.0320
		5	1.0021	6.57	24.0	
		6	1.0107	6.52	22.6	
		7	1.0028	6.52	23.2	
		8	1.0041	6.52	24.7	
		9	1.0134	6.51	31.3	
		10	1.0085	6.52	27.5	%Friable =

**Hardness readings in bold = tablets that capped during testing

TESTING PERFORMED BY A. FAVARA WEEK OF 03/13/06

N/A

To Page No. N/A

Recorded by: *John S. Deacony* Date: 03/21/06 Verified by: *QA* Date: 24 Oct 07

Project No. _____

Book No. _____

TITLE Physical Testing: Soma SR # 1038-24, 25, 26 03/21/06

From Page No. N/A

Physical Testing: Soma 700mg SR Tablets #1038-24, 25, & 26

Batch #	Target Wt. (mg)	Tablet #	Weight (g)	Thickness (mm)	Hardness (SCU)	Friability
1038-24	1000	1	1.0067	6.52	27.3	Initial Wt (g) = 10.0910 Final Wt (g) = 10.0500 %Friable =
		2	1.0081	6.55	24.7	
		3	1.0073	6.53	26.7	
		4	1.0097	6.56	26.3	
		5	1.0098	6.55	23.2	
		6	1.0046	6.54	25.9	
		7	1.0047	6.53	27.6	
		8	1.0051	6.53	24.8	
		9	0.9978	6.52	25.6	
		10	1.0206	6.55	27.5	
1038-25	1000	1	1.0078	6.71	26.7	Initial Wt (g) = 10.0490 Final Wt (g) = 10.0180 %Friable =
		2	1.0032	6.67	28.1	
		3	0.9992	6.68	28.3	
		4	0.9980	6.67	26.6	
		5	0.9954	6.67	28.3	
		6	1.0051	6.68	29.2	
		7	1.0036	6.68	22.4	
		8	1.0043	6.68	26.6	
		9	1.0061	6.70	27.9	
		10	1.0300	6.68	27.6	
1038-26	1000	1	1.0061	6.53	27.3	Initial Wt (g) = 10.0480 Final Wt (g) = 9.9980 %Friable =
		2	1.0016	6.50	27.2	
		3	0.9998	6.50	29.2	
		4	1.0003	6.50	26.3	
		5	1.0061	6.52	28.2	
		6	1.0020	6.50	26.6	
		7	1.0074	6.51	29.1	
		8	0.9928	6.47	25.6	
		9	0.9984	6.48	26.0	
		10	1.0036	6.52	25.1	

**Hardness readings in bold = tablets that capped during testing

00000000000000000000000000

NR 10
 XM 27.0 SC
 SREL 5.2 %
 SD 1.5 SC
 XMIN 25.1 SC
 XMAX 29.2 SC

10 25.1 SC
 9 26.0 SC
 8 25.6 SC
 7 29.1 SC
 ZERO
 6 26.6 SC
 ZERO
 5 28.2 SC
 ZERO
 4 26.3 SC
 ZERO
 3 29.2 SC
 ZERO
 2 27.2 SC
 1 27.3 SC

00000000000000000000000000



~~N/A~~

To Page No. N/A

Recorded by: John S. [Signature] Date: 03/21/06 Verified by: QVA Date: 24 Oct 07

Project No. _____

Book No. _____

TITLE ASTELIN 03-33 REGISTRATION BATCHES 03/21/06 99

From Page No. N/A

Azelastine Registration Batches - Decatur 03/2006

Batch # 03-33-03c
DOM: 03/13/06
Batch Size: 3300L

No problems occurred during the making of this batch. This batch was made in the 940 gallon processing tank, samples were taken according to E-Req 1112 sampling record (see below), the batch was filtered through a Millipore 0.45 micron PVDF 20" filter (code # CVHL-72P-P3), and transferred into the 940 gallon holding tank.

John S. [Signature] 03/21/06

E-Requisition 1112

SAMPLING RECORD

03-33-05s
03-33-04s

Batch Number: 03-33-03c

Bulk Solution

STEP	LOCATION	AMOUNT REQUIRED	PURPOSE	COLLECTED BY/DATE
PROCESSING Following Step 11	Top of Compounding Tank	2 x 300 mL	Azelastine and BAC Assay	N/A
	Bottom of Compounding Tank	2 x 300 mL	Azelastine and BAC Assay	
PROCESSING Following Step 18	Top of Holding Tank	2 x 300 mL	Azelastine and BAC Assay	
	Bottom of Holding Tank	2 x 300 mL	Azelastine and BAC Assay	
PROCESSING Final Bulk Solution Prior to Packaging		1 x 300 mL	Specific Gravity	
		1 x 500 mL	Spray Pattern Test Solution	
	Top of Holding Tank	approx. 220 mL	Microbial Limits	
		2 x 500 mL	Bulk Release Testing	
		2 x 500 mL (glass bottles)	Bulk Retain	

The only point of note is that during the filtration process it became clear that the filtration rate slowed considerably after approximately the 1-hour mark. At this time, I began charting the filtration rate (see pg. 107), and will continue to do so for the three remaining batches, for comparison.

John S. [Signature] 03/21/06

To Page No. 100

Recorded by: John S. [Signature] Date: 03/21/06 Verified by: QVA Date: 24 Oct 07

PTX0142-00009

Project No. _____

Book No. _____

TITLE ASTELIN 03-33 REGISTRATION BATCHES 03/21/06

From Page No. 99

Project No. _____

Book No. _____

TITLE ASTELIN 03-33 REGISTRATION BATCHES 03/21/06

101

03/21/06

Batch # 03-33-04s
DOM: 03/14/06
Batch Size: 330L

During the compounding of this batch, specifically steps # 8 through 11, it was observed that the 250 gallon processing tank began to shake (as it is equipped with shock-absorbing technology to counteract the force of the moving liquid inside) at a speed that caused the batch to "roll" or slosh around the inside of the tank, thus causing a loss of any true vortex. Without baffles to aid in the breaking up of the active and excipients, it was determined that the mixing speeds should be adjusted upwards until a suitable vortex could be obtained again.

Once this was done, no further problems arose in the manufacturing of the batch.

Samples were taken according to E-Req 1112 sampling record. The batch was filtered through a Millipore 0.45 micron PVDF 20" filter (code # CVHL-72P-P3), and transferred into a 375 gallon holding tank.

These new (greater) mixing speeds will be used for the two remaining 330L batches as well.

Batch # 03-33-05s
DOM: 03/15/06

No problems occurred during the making of this batch. This batch was made in the 250 gallon processing tank, samples were taken according to E-Req 1112 sampling record, the batch was filtered through a Millipore 0.45 micron PVDF 20" filter (code # CVHL-72P-P3), and transferred into a 375 gallon holding tank.

Batch # 03-33p-03s
DOM: 03/16/06

No problems occurred during the making of this batch. This batch was made in the 250 gallon processing tank, samples were taken according to E-Req 1113 sampling record (see below), the batch was filtered through a Millipore 0.45 micron PVDF 20" filter (code # CVHL-72P-P3), and transferred into a smaller holding tank.

ATTACHMENT D
SAMPLING RECORD

Batch Number: 03-33p-03s

Bulk Solution

ACTIVITY	SAMPLE LOCATION	AMOUNT REQUIRED	PURPOSE	COLLECTED BY/DATE
PROCESSING Final Bulk Solution	Top of holding tank	approx. 220 mL	Microbial Limits	* /
		2 x 500 mL	Release Testing	
		8 x 500 mL	Bulk Retain	

To Page No. 101

Recorded by: John S. Stewart Date: 03/21/06 Verified by: Em Bahwani Date: 10/25/07

Astelin 03-33 Registration Batches

Batch # 03-33-03c
D of M: 3/13/2006

Time	Batch Weight (kg)	Δ weight (kg)	
14:18	3350	n/a	after qs
15:20	214	3136	
15:21	200	14	
15:22	186	14	
15:23	174	12	
15:24	162	12	
15:25	150	12	
15:26	138	12	
15:27	126	12	
15:28	116	10	
15:29	104	12	
15:30	94	10	
15:31	84	10	
15:32	74	10	
15:33	64	10	
15:34	56	8	
15:35	48	8	
15:36	40	8	
15:37	30	10	
15:38	24	6	
15:39	14	10	
15:40	6	8	
15:41	4	2	
15:42	4	0	

Batch # 03-33-04s
D of M: 3/14/2006

Time	Batch Weight (kg)	Δ weight (kg)	
13:10:00	335	n/a	after qs
13:10:30	323	12	filter purge
13:11:00	313	10	
13:11:30	308	5	
13:12:00	302	8	
13:12:30	287	15	
13:13:00	265	22	
13:13:30	240	25	
13:14:00	213	27	
13:14:30	186	27	
13:15:00	159	27	
13:15:30	131	28	
13:16:00	103	28	
13:16:30	75	28	
13:17:00	47	28	
13:17:30	18	29	
13:18:00	0	18	

Batch # 03-33-05s
D of M: 3/15/2006

Time	Batch Weight (kg)	Δ weight (kg)	
	335	n/a	after qs
12:35:00	328	7	after sampling
12:35:30	308	20	after filter purge
12:36:00	270	38	
12:36:30	224	46	
12:37:00	178	46	
12:37:30	133	45	
12:38:00	88	45	
12:38:30	43	45	
12:39:00	2	41	
12:39:30	2	0	
12:40:00	0	2	

Batch # 03-33p-03s
D of M: 3/16/2006

Time	Batch Weight (kg)	Δ weight (kg)	
	335	n/a	after qs
12:42:00	306	29	after filter purge
12:42:30	292	14	
12:43:00	248	44	
12:43:30	203	45	
12:44:00	158	45	
12:44:30	114	44	
12:45:00	71	43	
12:45:30	37	34	
12:49:00*	34	3	
12:49:30	3	31	
12:50:00	3	0	

* FILTRATION WAS PAUSED AS THE HOLDING TANK WAS NEARING CAPACITY. WE ESTIMATED THAT THERE WAS ENOUGH ROOM & RESUMED FILTRATION.

John S. Stewart 03/21/06
N/A

03/21/06
JS
101

To Page No. 101

Recorded by: John S. Stewart Date: 03/21/06 Verified by: Em Bahwani Date: 10/25/07

Project No. _____

Book No. _____

TITLE TUSSI-ORGANIDIN NR LIQUID

04/04/06

From Page No. N/A

Product: TUSSI-ORGANIDIN NR LIQUID

Form #: 1044-102P

Re: Codeine placebo of Tussi-Organidin NR for analytical testing.

BATCH SIZE: 2.000 L theoretical SG = 1.125

#	Ingredient	Code #	Lot #	mg per 5mL	theor. grams per batch (g)	actual grams per batch (g)
1	Guaifenesin USP	4055	000002356	100.0	40.00	40.014
3	Sodium Benzoate NF	0168	2L 2508	5.0	2.00	2.007
4	Saccharin Sodium USP	7260	000001994	2.0	0.80	0.806
5	Citric Acid Anhydrous USP	1134	000001292	4.6	1.84	1.848
6	Propylene Glycol USP	6257	0D0508	500.0	200.00	200.2
7	Glycerin USP	5100	000002135	2160.0	863.50	863.5
8	Sorbitol Solution 70% USP	8135	000001230	1000.0	400.00	400.0
9	Dye: FD&C Red #40	3073	000002602	1.0	0.40	0.402
10	Flavor: Raspberry	3425	000002416	18.7	7.48	7.487
11	Water, purified	8542	Q.S. (to make 2249.4g)			

TOTAL: 2249.4 g

Theoretical specific gravity is 1.125 based on Processing Batch Record

Batch produced in Room #133

Manufacturing Procedure:

- In a stainless steel beaker (FINAL CONTAINER 724.8 g tare), add approximately 600g of Purified water.
- With mixing, add the Red #40 (✓) to the beaker and mix until uniform.
- Lower the mixing speed and slowly add the Propylene Glycol (✓) and mix until uniform.
- Raise the mixing speed again to obtain a sufficient vortex and add the Guaifenesin (✓), mixing until the solution is complete. Mixing Time: 5 min.
- Lower the mixing speed and add the Glycerin (✓) and Sorbitol Solution (✓). Mix until uniform. Mixing Time: 5 min.
- Raise the mixing speed and slowly add the following, mixing after each until the solution is complete:
 - Sodium Benzoate (✓)
 - Saccharin Sodium (✓)
 - Citric Acid (✓)
 - Raspberry Flavor (✓)
- Adjust the batch weight to 2249.4g with Purified water.

Before adjustment:	After Adjustment:
Gross: <u>3097.6</u>	Gross: <u>N/A</u>
Tare: <u>724.8</u>	Tare: <u>X</u>
Net: <u>2372.8</u>	Net: <u>X</u>

Amount of water needed _____ g to bring final net weight to 2249.4g

- Mix for 10 minutes or until uniform.

Mixing Time: _____, Final pH: _____ (date _____), Final Viscosity: _____ (date _____)

Viscosity Testing - spindle # _____, 60 RPM for 1 minute @ 25°C

JSD 04/04/06

BATCH HAS BEEN OVER-DILUTED WITH RINSE WATER DURING MANUFACTURING. WILL BE REPEATED

To Page No. N/A

Recorded by: John S. [Signature]

Date: 04/04/06

Verified by: Enl Bahwani

Date: 10/25/07

Project No. _____

Book No. _____

TITLE TUSSI-ORGANIDIN NR LIQUID

04/05/06

103

From Page No. N/A

Product: TUSSI-ORGANIDIN NR LIQUID

Form #: 1044-103P

Re: Codeine placebo of Tussi-Organidin NR for analytical testing.

BATCH SIZE: 2.000 L theoretical SG = 1.125

#	Ingredient	Code #	Lot #	mg per 5mL	theor. grams per batch (g)	actual grams per batch (g)
1	Guaifenesin, USP	4055	000002356	100.0	40.00	40.009
3	Sodium Benzoate, NF	0168	2L 2508	5.0	2.00	2.007
4	Saccharin Sodium, USP	7260	000001994	2.0	0.80	0.802
5	Citric Acid, USP, Anhydrous	1134	000001292	4.6	1.84	1.842
6	Propylene Glycol, USP	6257	0D0508	500.0	200.00	200.0
7	Glycerin, USP	5100	000002135	2160.0	863.50	863.5
8	Sorbitol Solution, USP, 70%	8135	000001230	1000.0	400.00	400.0
9	Dye: FD&C Red #40	3073	000002602	1.0	0.40	0.405
10	Flavor: Raspberry	3425	000002416	18.7	7.48	7.484
11	Water, purified	8542	Q.S. (to make 2249.4g)			

TOTAL: 2249.4 g

Theoretical specific gravity is 1.125 based on Processing Batch Record

Batch produced in Room #133

Manufacturing Procedure:

- In a stainless steel beaker (FINAL CONTAINER 598.5 g tare), add approximately 600g of Purified water.
- With mixing, add the Red #40 (✓) to the beaker and mix until uniform.
- Lower the mixing speed and slowly add the Propylene Glycol (✓) and mix until uniform.
- Raise the mixing speed again to obtain a sufficient vortex and add the Guaifenesin (✓), mixing until the solution is complete. Mixing Time: 5 min.
- Lower the mixing speed and add the Glycerin (✓) and Sorbitol Solution (✓). Mix until uniform. Mixing Time: 10 min.
- Raise the mixing speed and slowly add the following, mixing after each until the solution is complete:
 - Sodium Benzoate (✓)
 - Saccharin Sodium (✓)
 - Citric Acid (✓)
 - Raspberry Flavor (✓)
- Adjust the batch weight to 2249.4g with Purified water.

Before adjustment:	After Adjustment:
Gross: <u>2720.5</u>	Gross: <u>2847.9</u>
Tare: <u>598.5</u>	Tare: <u>598.5</u>
Net: <u>2122.0</u>	Net: <u>2249.4</u>

Amount of water needed 127.4 g to bring final net weight to 2249.4g

- Mix for 10 minutes or until uniform.

Mixing Time: 20 min., Final pH: 4.25 (date 04/05/06), Final Viscosity: 18.0 sp (date 04/05/06)

Viscosity Testing - spindle # 2, 60 RPM for 1 minute @ 25°C

John S. [Signature] 04/05/06

To Page No. N/A

Recorded by: John S. [Signature]

Date: 04/05/06

Verified by: Enl Bahwani

Date: 10/25/07

Project No. _____

Book No. _____

TITLE ASTELIN 03-33 REG. BATCHES

03/24/06

From Page No. 104

Astelin 03-33 Registration Batches

John S. Balwani 03/24/06

Batch # 03-33-03c
DOM: 3/13/2006

Time	Batch Weight (kg)	Amount Filtered (kg)	Δ weight (kg)	Rate (kg/min)
14:18	3350	n/a	n/a	n/a
15:20	214	3136	3136	50.6
15:21	200	3150	14	14.0
15:22	186	3164	14	14.0
15:23	174	3176	12	12.0
15:24	162	3188	12	12.0
15:25	150	3200	12	12.0
15:26	138	3212	12	12.0
15:27	126	3224	12	12.0
15:28	116	3234	10	10.0
15:29	104	3246	12	12.0
15:30	94	3256	10	10.0
15:31	84	3266	10	10.0
15:32	74	3276	10	10.0
15:33	64	3286	10	10.0
15:34	56	3294	8	8.0
15:35	48	3302	8	8.0
15:36	40	3310	8	8.0
15:37	30	3320	10	10.0
15:38	24	3326	6	6.0
15:39	14	3336	10	10.0
15:40	6	3344	8	8.0
15:41	4	3346	2	2.0
15:42	4	3346	0	0.0

Batch # 03-33-04s
DOM: 3/14/2006

Time	Batch Weight (kg)	Amount Filtered (kg)	Δ weight (kg)	Rate (kg/min)
13:10:00	335	n/a	n/a	n/a
13:10:30	323	12	12	n/a
13:11:00	313	22	10	22.0
13:11:30	308	27	5	15.0
13:12:00	302	33	6	11.0
13:12:30	287	48	15	21.0
13:13:00	265	70	22	37.0
13:13:30	240	95	25	47.0
13:14:00	213	122	27	52.0
13:14:30	186	149	27	54.0
13:15:00	159	176	27	54.0
13:15:30	131	204	28	55.0
13:16:00	103	232	28	56.0
13:16:30	75	260	28	56.0
13:17:00	47	288	28	56.0
13:17:30	18	317	29	57.0
13:18:00	0	335	18	47.0

N/A

To Page No. 105

Recorded by: *John S. Balwani*

Date: 03/24/06

Verified by: *Em Balwani*

Date: 10/25/07

TITLE ASTELIN 03-33 REGISTRATION BATCHES

Project No. _____

Book No. _____

03/24/06 105

From Page No. 104

John S. Balwani 03/24/06

Batch # 03-33-05s
DOM: 3/15/2006

Time	Batch Weight (kg)	Amount Filtered (kg)	Δ weight (kg)	Rate (kg/min)
	335	n/a	n/a	n/a
12:35:00	328	n/a	7	n/a
12:35:30	308	n/a	20	n/a
12:36:00	270	65	38	58.0
12:36:30	224	111	46	84.0
12:37:00	178	157	46	92.0
12:37:30	133	202	45	91.0
12:38:00	88	247	45	90.0
12:38:30	43	292	45	90.0
12:39:00	2	333	41	86.0
12:39:30	2	333	0	41.0
12:40:00	0	335	2	2.0

Batch # 03-33p-03s
DOM: 3/16/2006

Time	Batch Weight (kg)	Amount Filtered (kg)	Δ weight (kg)	Rate (kg/min)
	335	n/a	n/a	n/a
12:42:00	306	n/a	29	n/a
12:42:30	292	43	14	43.0
12:43:00	248	87	44	58.0
12:43:30	203	132	45	89.0
12:44:00	158	177	45	90.0
12:44:30	114	221	44	89.0
12:45:00	71	264	43	87.0
12:45:30	37	298	34	77.0
12:49:00	34	301	3	37.0
12:49:30	3	332	31	34.0
12:50:00	3	332	0	31.0

N/A

To Page No. 106

Recorded by: *John S. Balwani*

Date: 03/24/06


Verified by: *Em Balwani*

Date: 10/25/07

PTX0142-00012

Filtration Study: Astelin Registration Batches
(Decatur - March 2006)

John S. Wilcox

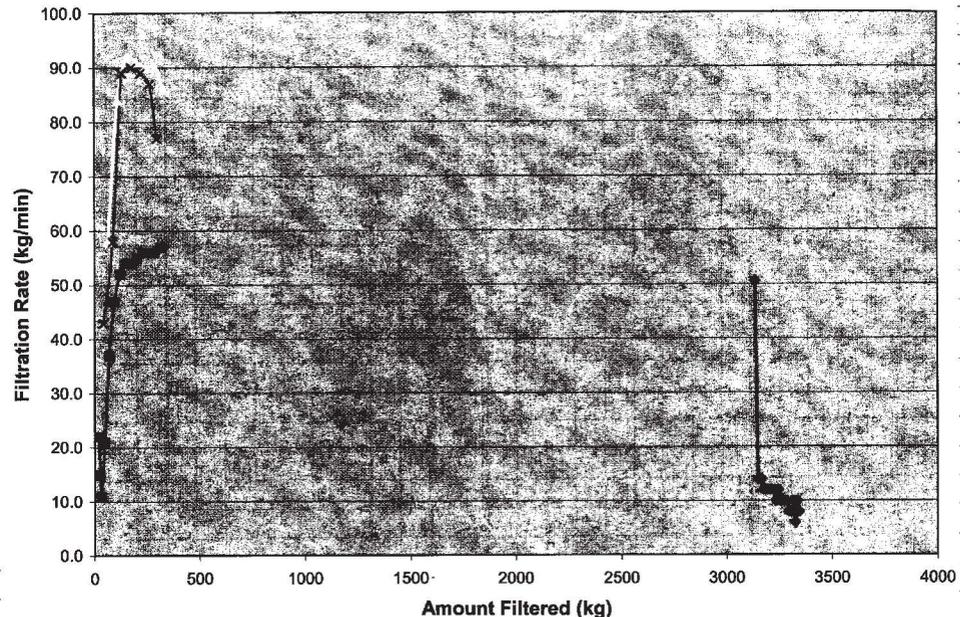

◆ 03-33-03c ■ 03-33-04s 03-33-05s ✕ 03-33-p-03s

N/A

To Page No. N/A

Recorded by: *John S. Wilcox* Date: 03/24/06 Verified by: *Em Bahwani* Date: 10/25/07

From Page No. N/A

AT THE REQUEST OF C. YAMAC, 20 BOTTLES OF ASTELIN REFORMULATED #1027-17 WERE PLACED ON FREEZE/THAW STABILITY.

TEN BOTTLES WERE NUMBERED & PLACED UPRIGHT
TEN BOTTLES WERE NUMBERED & PLACED IN THE PRONE POSITION.

THEY WERE PLACED IN THE FREEZER IN ROOM #129 AND WILL UNDERGO FREEZE/THAW WITH EXTREME DIFFERENCES (FROM FREEZER TO 40°C OVEN).

WED. 03/29/06: 20 SAMPLES MOVED FROM FREEZER IN ROOM #129 TO 40°C/AMBIENT OVEN IN STABILITY CHAMBER ROOM = THAW #1.

FRI. 03/31/06: 20 SAMPLES, WHICH HAVE THAWED, ARE REMOVED FROM THE 40°C OVEN AND RETURNED TO THE FREEZER IN ROOM #129. NO CHANGES TO THE OUTSIDE OF THE CONTAINERS OBSERVED.

MON. 04/03/06: 20 SAMPLES ARE MOVED FROM THE FREEZER IN Rm #129, BACK TO THE 40°C/AMBIENT OVEN IN THE STABILITY CHAMBER ROOM FOR THAW #2.

WED. 04/05/06: 20 SAMPLES ARE MOVED FROM 40°C/AMBIENT OVEN, TO THE FREEZER IN ROOM #129 FOR FREEZE #3. NO CHANGES TO THE OUTSIDE OF THE CONTAINERS WAS OBSERVED.

FRI. 04/07/06: 20 SAMPLES REMOVED FROM FREEZER & PLACED IN 40°C/AMBIENT. NO CHANGES OBSERVED UPON FREEZE #3.

MON. 04/10/06: 20 SAMPLES REMOVED FROM 40°C/AMBIENT AFTER THAW #3. NO CHANGES OBSERVED TO THE OUTSIDE OF THE CONTAINER.

BOTTLES WILL BE PLACED IN STABILITY LAB FOR RETAIN/Possible FUTURE TESTING. JSD 04/11/06

To Page No. N/A

Recorded by: *John S. Wilcox* Date: 03/27/06 Verified by: *Em Bahwani* Date: 10/25/07

Project No. _____

Book No. _____

TITLE ASTELIN NASAL SPRAY
FUSST-ORGANIZATION NR JSD 04/05/06 04/07/06
FORM. 03-33

From Page No. 111

Product: Astelin® Nasal Spray

pg. 1044-110

Re: Lab batch to match Investigational Formulation 03-33, for Analytical Testing.

Ingredient	Code #	Lot #	%w/v	Theo. grams per batch	Actual grams per batch
Azelastine HCl	1105	000002047	0.100%	2.000	2.011
HPMC 2910, USP	3340	000002822	0.100%	2.000	2.002
Sorbitol Solution, USP, 70%	8135	000001230	6.446%	128.920	128.9
Sucralose, NF	7474	04073	0.150%	3.000	3.008
Edetate Disodium, USP	1349	000002353	0.050%	0.500	0.502
Sodium Citrate, USP, Dihydrate	0734	2C1804	0.068%	1.360	1.368
Benzalkonium Chloride, NF, 50% Soln.	2202	1C1304	0.025%	0.500	0.502
Water, Purified or Deionized			QS to 100%	QS to 2030.0g	

Manufacturing Procedure

- To a suitably sized container (tare weight 726.7 g), obtain 740g of deionized water. Begin mixing and heat to 80°C (175°F). Record the temperature. While mixing, slowly add the HPMC (✓) to the vortex of the mixing water. Mix for 20 minutes or until the mixture is lump-free. Record the mixing time. Do NOT let the mixture cool to below 60°C (140°F) while mixing. Temp: 80 °C Mixing Time: 20 min.
- Add approx. 1000g of deionized water. Cool to 32°C (90°F). Record the temperature. Temp: 28 °C
- Add the Azelastine HCl (✓), mixing until completely dissolved. Record the mixing time and temperature. Temp: 28 °C Mixing Time: 10 min.
- Add the following ingredients, one at a time, to the mixture: sucralose (✓), sorbitol solution (✓), EDTA (✓), sodium citrate (✓). Mix after each until solution is clear and record the overall mixing time. Mixing Time: 20 min.
- Add the Benzalkonium Chloride (✓) to the batch and mix for 10 minutes. Mixing Time: 10 min.
- Adjust batch weight to (2030.0g) with deionized water. Mix for 30 minutes or until uniform.

Before adjustment:	After adjustment:
Gross: <u>2673.0 g</u>	Gross: <u>2756.7 g</u>
Tare: <u>726.7 g</u>	Tare: <u>726.7 g</u>
Net: <u>1946.3 g</u>	Net: <u>2030.0 g</u>

Amount of Water added: 83.7 g to bring final net weight to 2030.0g.

Mixing Time: 30 min. Temperature at end of mixing: 22 °C Description: CLEAR, VERY

- Filter finished product through a 0.22 µm Millipore Stainless Steel High Pressure Filter Holder. LITTLE FOAM
Filter Info: MILLIPORE 0.22 µm GVWP14250 LOT# ROAN 29130
Discard the first 50mL to 100mL and retain the rest of the product for testing and packaging.

- Measure the pH, viscosity, and specific gravity (viscosity testing – spindle #2, 60rpm for 1 minute at 25°C).

Final pH: 6.36 (date 04/10/06). Final viscosity: 5.0 cps (date 04/10/06).

Specific Grav.: _____ (date _____) [pyc.: _____ g / pyc+water = _____ g / pyc+sample = _____ g]

To Page No. 111

Recorded by: John S. [Signature]

Date: 04/07/06

Verified by: Enl Bahwani

Date: 10/25/07

Project No. _____

Book No. _____

TITLE ASTELIN NASAL SPRAY
FUSST-ORGANIZATION NR JSD 04/05/06 04/07/06
FORM. 03-33p

111
04/07/06

From Page No. 111

Product: Astelin® Nasal Spray

pg. 1044-111

Re: Lab batch to match Investigational Placebo 03-33p, for Analytical Testing.

Ingredient	Code #	Lot #	%w/v	Theo. grams per batch	Actual grams per batch
Azelastine HCl	1105		0.000%	0.000	
HPMC 2910, USP	3340	000002822	0.100%	2.000	2.006
Sorbitol Solution, USP, 70%	8135	000001230	6.446%	128.920	128.9
Sucralose, NF	7474	04073	0.150%	3.000	3.009
Edetate Disodium, USP	1349	000002353	0.050%	0.500	0.504
Sodium Citrate, USP, Dihydrate	0734	2C1804	0.068%	1.360	1.361
Benzalkonium Chloride, NF, 50% Soln.	2202	1C1304	0.025%	0.500	0.501
Water, Purified or Deionized			QS to 100%	QS to 2028.0g	

Manufacturing Procedure

- To a suitably sized container (tare weight 724.9 g), obtain 740g of deionized water. Begin mixing and heat to 80°C (175°F). Record the temperature. While mixing, slowly add the HPMC (✓) to the vortex of the mixing water. Mix for 20 minutes or until the mixture is lump-free. Record the mixing time. Do NOT let the mixture cool to below 60°C (140°F) while mixing. Temp: 78 °C Mixing Time: 20 min.
- Add approx. 1000g of deionized water. Cool to 32°C (90°F). Record the temperature. Temp: 32 °C
- ~~Add the Azelastine HCl (✓), mixing until completely dissolved. Record the mixing time and temperature. Temp: _____ °C Mixing Time: _____ min.~~
PLACEBO Temp: _____ °C Mixing Time: _____ min.
- Add the following ingredients, one at a time, to the mixture: sucralose (✓), sorbitol solution (✓), EDTA (✓), sodium citrate (✓). Mix after each until solution is clear and record the overall mixing time. Mixing Time: 20 min.
- Add the Benzalkonium Chloride (✓) to the batch and mix for 10 minutes. Mixing Time: 15 min.
- Adjust batch weight to (2028.0g) with deionized water. Mix for 30 minutes or until uniform.

Before adjustment:	After adjustment:
Gross: <u>2623.0 g</u>	Gross: <u>2702.9 g</u>
Tare: <u>724.9 g</u>	Tare: <u>724.9 g</u>
Net: <u>1898.1 g</u>	Net: <u>2028.0 g</u>

Amount of Water added: 129.9 g to bring final net weight to 2028.0g.

Mixing Time: 30 min. Temperature at end of mixing: 28 °C Description: CLEAR, SOME

- Filter finished product through a 0.22 µm Millipore Stainless Steel High Pressure Filter Holder. FOAM
Filter Info: MILLIPORE 0.22 µm GVWP14250
Discard the first 50mL to 100mL and retain the rest of the product for testing and packaging. Lot # ROAN 29130

- Measure the pH, viscosity, and specific gravity (viscosity testing – spindle #2, 60rpm for 1 minute at 25°C).

Final pH: 6.27 (date 04/10/06). Final viscosity: 5.0 cps (date 04/10/06).

Specific Grav.: _____ (date _____) [pyc.: _____ g / pyc+water = _____ g / pyc+sample = _____ g]

To Page No. 111

Recorded by: John S. [Signature]

Date: 04/07/06

Verified by: Enl Bahwani

Date: 10/25/07

Project No. _____

Book No. _____

TITLE OSMOLARITY TESTING: OPTIVAR 04/14/06
LOT # 02306A

From Page No. 111

OPTIVAR (AZELASTINE HCL EYE DROPS 0.05%) Lot #02306A was received from CARDINAL, from a 250 L BATCH MANUFACTURED ON 04/04/06. IT WILL UNDERGO OSMOLARITY TESTING, AN OSMOLARITY ADJUSTMENT USING SORBITOL SOLUTION, AND REPEAT OSMOLARITY TESTING TO OBSERVE ANY CHANGE.

THE OSMETTE-A AUTOMATIC OSMOMETER (MODEL #5002, SER #11110073) WAS CALIBRATED AS FOLLOWS:

- TURN TOGGLE TO "OPERATE"
- ALLOW MACHINE TO WARM UP FOR ≈ 15 MIN.
- ADD COOLANT UNTIL FULL
- PREPARE MULTIPLE 500 mOsm/kg & 100 mOsm/kg STD. SAMPLES.
- SET RANGE TO 0-2
OPERATE TO SPECIAL
- PLACE 500 STD INTO MACHINE, ADJUST READINGS WITH "I" KNOB.
- SWITCH TO OPERATE AND READ A NEW, CLEAN 500 STD.
RESULT = 504, REPEAT: RESULT = 501
- REPEAT WITH 100 mOsm/kg STD, ADJUSTING WITH "II" KNOB.
RESULT = 97, 101
- TEST OPTIVAR SAMPLE (CARDINAL RESULT ≈ 250 mOsm/kg)
RESULT #1 = 242 mOsm/kg
#2 = 241 mOsm/kg

OSMOMETER CALIBRATED... SEE PG. 113 FOR DIRECTIONS FOR OSMOLARITY ADJUSTMENT OF SAMPLE (LOT # 02306A) AND RESULTS.

~~N/A~~

To Page No. 113

Recorded by: <i>John S. [Signature]</i>	Date: <u>04/14/06</u>	Verified by: <i>Enl Balwani</i>	Date: <u>10/25/07</u>
---	-----------------------	---------------------------------	-----------------------

Project No. _____

Book No. _____

TITLE OSMOLARITY TESTING: OPTIVAR 04/14/06

From Page No. 112

ADJUSTMENT #1

Adjustment of Osmolarity of Sorbitol Solution
Azelastine Hydrochloride Eye Drops 0.05% Lot 02306A, 250 liter batch made at Cardinal. Date 04/04/2006

Amount of Sorbitol Solution Required for the batch = 16666.5g
 Amount of Sorbitol Solution Added in the batch = 16000.5g
 Difference = 666.5g/250 liter
 = 666.5g/253.75 kg
 = 2.6266g/kg
 = 0.0026266g/g

To 10g of the product add 0.02627g of Sorbitol solution.
Measure Osmolarity.

John S. [Signature] 04/14/06

LOT # 02306A = 10.005g
SORBITOL SOL 70% = 0.0274g (Lot # 0000001230)

OSMOLARITY RESULTS:

TEST # 1 = 251 mOsm/kg
TEST # 2 = 251 mOsm/kg

ADJUSTMENT #2

WEIGHED 5x SORBITOL SOLUTION 70% & QS TO 50g TOTAL. (THEO. SOLB. = 0.13135g)

SORBITOL SOLUTION, 70% (Lot # 000001230) = 0.1390g (ACTUAL)
OPTIVAR (Lot # 02306A) QS TO 50g = 50.0048g (ACTUAL)

RETEST 500 mOsm/kg STD = 502
RETEST 100 mOsm/kg STD = 103
ADJUSTED OPTIVAR #2a = 251
ADJUSTED OPTIVAR #2b = 251 mOsm/kg

To Page No. N/A

Recorded by: <i>John S. [Signature]</i>	Date: <u>04/14/06</u>	Verified by: <i>Enl Balwani</i>	Date: <u>10/25/07</u>
---	-----------------------	---------------------------------	-----------------------

Project No. _____

Book No. _____

TITLE ASTELIN + STEROID EXPERIMENT 04/24

From Page No. 114

FLUTICASONA PROPIONATE NASAL SPRAY (LOT #C184889)

17ml AZELASTINE HCL NS COMMERCIAL BOTTLE + PUMP

= 10.9713g

BOTTLE w/o PUMP

= 4.5717g

COMMERCIAL BOTTLE + FLUTICASONA N.S.

= 19.1669g

AMOUNT OF FLUTICASONA N.S. ADDED

= 14.5952g

AMOUNT OF AZELASTINE HCL NEEDED FOR 0.1% w/v

= 0.0146g

AMOUNT OF AZELASTINE HCL WEIGHED:

= 0.0158g (0.0009g LEFT ON WEIGHBOAT ∴ AMOUNT ADDED = 0.0149g)

BOTTLE IS THEN SONICATED FOR 15 MINUTES, AND WILL BE OBSERVED.

FINAL SAMPLE APPEARS TO BE UNIFORM, WITH NO PRECIP.

N/A

To Page No. 116

Recorded by: *[Signature]*

Date: 04/24/06

Verified by: *[Signature]*

Date: 10/25/07

Project No. _____

Book No. _____

04/25/06 115

STABILITY TESTING: Soma SR #1048-10, 12, 14
66, 82

Stability Testing: Soma 700mg SR Tablets #1048-10, 12, 14 @ 2 months (40°C/75%RH)

Lot #	Target Wt. (mg)	Tablet #	Weight (g)	Thickness (mm)	Hardness (SCU)	Friability
1048-10 with desiccant	1000	1	0.9827	6.62	18.0	Initial Wt (g) = 10.0080
		2	1.0075	6.99	21.9	Final Wt (g) = 9.9590 %Friable = 21.6
		3	1.0009	6.95	19.3	
		4	1.0055	6.97	21.2	
		5	0.9899	6.67	21.6	

Lot #	Target Wt. (mg)	Tablet #	Weight (g)	Thickness (mm)	Hardness (SCU)	Friability
1048-12 with desiccant	1000	1	0.9964	6.80	19.1	Initial Wt (g) = 9.9660
		2	0.9855	6.77	17.2	Final Wt (g) = 9.9160 %Friable = 16.8
		3	1.0117	6.85	22.9	
		4	0.9944	6.77	18.4	
		5	0.9849	6.74	16.8	

Lot #	Target Wt. (mg)	Tablet #	Weight (g)	Thickness (mm)	Hardness (SCU)	Friability
1048-14 without desiccant	1000	1	0.9965	7.01	18.0	Initial Wt (g) = 9.9870
		2	0.9975	7.02	17.7	Final Wt (g) = 9.9380 %Friable = 17.8
		3	1.0021	7.02	16.6	
		4	1.0008	6.99	18.1	
		5	1.0072	7.06	17.8	

Stability Testing: Soma 700mg SR Tablets #1048-66 & 82 @ 2 months (40°C/75%RH)

Lot #	Target Wt. (mg)	Tablet #	Weight (g)	Thickness (mm)	Hardness (SCU)	Friability
1048-66 with desiccant	950	1	0.9436	6.45	25.1	Initial Wt (g) = 9.4710
		2	0.9330	6.45	24.7	Final Wt (g) = 9.4400 %Friable = 27.5
		3	0.9489	6.45	23.8	
		4	0.9484	6.45	29.7	
		5	0.9585	6.48	27.5	

Lot #	Target Wt. (mg)	Tablet #	Weight (g)	Thickness (mm)	Hardness (SCU)	Friability
1048-66 without desiccant	950	1	0.9558	6.49	26.6	Initial Wt (g) = 9.3760
		2	0.9599	6.54	27.2	Final Wt (g) = 9.3380 %Friable = 26.3
		3	0.9482	6.48	24.8	
		4	0.9421	6.45	25.0	
		5	0.9641	6.63	26.3	

Lot #	Target Wt. (mg)	Tablet #	Weight (g)	Thickness (mm)	Hardness (SCU)	Friability
1048-82 with desiccant	1000	1	0.9999	6.75	26.3	Initial Wt (g) = 10.0050
		2	0.9958	6.68	25.9	Final Wt (g) = 9.9780 %Friable = 26.7
		3	1.0065	6.75	26.4	
		4	0.9985	6.73	27.3	
		5	0.9987	6.71	26.7	

Lot #	Target Wt. (mg)	Tablet #	Weight (g)	Thickness (mm)	Hardness (SCU)	Friability
1048-82 without desiccant	1000	1	1.0024	6.73	26.2	Initial Wt (g) = 9.9390
		2	0.9912	6.67	29.1	Final Wt (g) = 9.9180 %Friable = 26.7
		3	0.9901	6.67	26.7	
		4	0.9919	6.67	27.3	
		5	0.9957	6.69	26.7	

not enough tablets remaining to perform hardness testing.

10/25/07
[Signature]

To Page No. 114

Recorded by: *[Signature]*

Date: 04/25/06

Verified by: *[Signature]*

Date: 24 Oct 07

FUTICHSONE PROPIONATE NASAL SPRAY, 50 mg
 Lot # C184889 Ex. P. Rev. 2007

13.0834g = CENTRIFUGE TUBE THREE WT.
 14.5465g = FUTICHSONE N.S.

SAMPLE WAS PLACED IN CENTRIFUGE FOR 30 MIN. OPERATIONS
 WAS STILL CLOUDY, ∴ RETURNED TO CENTRIFUGE FOR AN
 ADDITIONAL 30 MINUTES.

THE SUPERNATANT WAS STILL CLOUDY. IT WAS FOUNDED OFF INTO
 A SCINTILLATION VIAL. THE SAMPLE WAS OPAQUE WHITE.

A MICROSCOPY SLIDE WAS PREPARED WITH A DROP OF THE SAMPLE
 AND VIEWED @ 40x MAGNIFICATION. THE SLIDE APPEARED TO
 BE FREE OF ANY SOLIDS, HOWEVER THE SAMPLE WAS PREPARED
 TO BE FILTERED.

USING A 10ml SPRING AND A 5ml STERILE ACRODISC, THE
 SAMPLE WAS FILTERED AND PLACED IN A NEW SCINTILLATION VIAL.

THE SAMPLE WAS NOT CLEAR (STILL OPAQUE AND WHITE) ∴ IT WAS
 PASSED THROUGH A 2µm STERILE ACRODISC. THE RESULTING
 SAMPLE WAS STILL OPAQUE ∴ IT WAS NEXT PASSED THROUGH A
 1µm FILTER. THE RESULTING SAMPLE WAS STILL OPAQUE ∴ IT
 WAS DECIDED THAT IT WOULD NEXT BE PASSED THROUGH A 0.45µm
 STERILE ACRODISC.

ATTEMPTING TO PASS THE SAMPLE THROUGH THE 0.45µm FILTER
 WAS EXTREMELY DIFFICULT. AFTER ~3ml OF THE SAMPLE HAD
 PASSED THROUGH (AND WAS SEEMINGLY CLEAR NOW) THE PRESSURE
 NECESSARY TO FORCE THE SAMPLE THROUGH THE FILTER WAS
 MORE THAN THE FILTER COULD WITHSTAND. THE ACRODISC GAVE
 WAY AROUND ITS CIRCUMFERENCE AND THE REMAINING SAMPLE
 WAS SPILLED AND LOST.
 THIS PROCESS WILL BE ATTEMPTED AGAIN TOMORROW.

Recorded by: [Signature] Date: 04/25/06 Verified by: Kim Babarum Date: 10/25/07 To Page No. 117

15.3786g = SCINTILLATION VIAL + STIR ROD

23.5168g = VIAL + STIRROD + PRODUCT *

8.1382g = PRODUCT (FUTICHSONE N.S. Lot # C184889)

∴ 0.0081g = AERUSIVE HCL (Lot # 0000002047) TO BE ADDED = 0.1%

∴ 23.5249g = TARGET GROSS WT.

23.52492g = ACTUAL WT. (gross)

* THE SAMPLE WAS PREPARED IN A SIMILAR FASHION TO THE SAMPLE FROM PAGES
 FUTICHSONE PROPIONATE NASAL SPRAY (Lot # C184889) SAMPLE WAS BLENDED
 SHAKE WELL, TO MIX AND SUSPEND THE SAMPLE PROPERLY. THE SPRAY BOTTLE
 WAS CUT OFF THE AMBER GLASS BOTTLE AND Poured INTO A 50ml CENTRIFUGE
 VIAL. THE SAMPLE WAS SHUN FOR ONE HOUR, THE SUPERNATANT POWDER
 OFF, THE SOLIDS WERE DISCARDED
 THE SUPERNATANT SAMPLE WAS INTRODUCED INTO A 10ml SPRING AND
 PASSED THROUGH AN UNCONNECTED SERIES OF 5µm, 2µm, AND 1µm
 STERILE ACRODISC (P.VDF) FILTERS. THE RESULTING SAMPLE WAS THEN
 SLOWLY AND DIFFICULTLY PASSED THROUGH A 0.45µm ACRODISC.

THIS SAMPLE WAS TRANSFERRED INTO A PRE-WEIGHED SCINTILLATION VIAL
 CONTAINING A METRIC STIR ROD, A 0.1% w/v SOLUTION OF AERUSIVE
 HCL WAS ADDED, USING AN ANALYTICAL BALANCE, AND THE SAMPLE WAS PLACED
 ON A STIRRED/HORIZONTAL AND MIXED OVERNIGHT.

THE FOLLOWING MORNING THE SAMPLE WAS OBSERVED FOR CLARITY/CLUMPING
 AND IT WAS FOUND THAT THERE WAS STILL A LARGE PORTION OF THE
 AERUSIVE HCL THAT HAD NOT GONE INTO SOLUTION.

THIS SAMPLE WAS COVERED (CAPPED) AND LABELLED FOR POSSIBLE FURTHER
 EXPERIMENTATION.

Recorded by: [Signature] Date: 04/26/06 Verified by: Kim Babarum Date: 10/25/07 To Page No. 118

THE REMAINING SAMPLE FROM YESTERDAY'S EXPERIMENT WAS PLACED ON A SPINNING/HOTPLATE. A THERMOMETER WAS PLACED IN THE SAMPLE.

INITIAL TEMP. = 19°C
 TARGET TEMP. = 60°C OR OBSERVED CAPACITY OF THE SAMPLE.

SAMPLE WAS HEATED, WITH STIRRING, TO $\approx 60^\circ\text{C}$ (THERMOMETER = 69°C) AND HEAT STOPPED FOR 15 MIN. AFTER 15 MIN. THE STIRRING WAS STOPPED AND THE SAMPLE WAS REMOVED FROM THE HEAT.

THE SAMPLE REMAINED TRANSLUCENT AND IT IS ASSUMED THAT IF ALLOWED TO SIT LONG ENOUGH, THE UNDISSOLVED AZEOSTRINE WILL BECOME VISIBLE AGAIN.

AFTER 30 MIN, NO SOLID WAS VISIBLE. THE NEW WHITE OPAQUE SAMPLE WAS PLACED BACK ON THE STIRRER/HOTPLATE.

INITIAL TEMP. = 21°C
 TARGET TEMP. = 60°C FOR ≈ 15 MORE MIN, OR CLARITY OF SAMPLE.

AFTER ANOTHER 15 MINUTES OF HEATING, THE SAMPLE DID NOT CHANGE. IT REMAINS WHITE AND OPAQUE. IT WILL BE ALLOWED TO SETTLE AND FULLY COOL & BE OBSERVED FOR PRECIP. 04/27/06

THE FOLLOWING MORNING (04/28/06) IT IS OBSERVED THAT SOME OF THE SOLID HAS, IN FACT, PRECIPITATED OUT OF THE 1044-117 SAMPLE.

IT IS DEEMED NECESSARY TO SUBMIT 1044-117 TO ANALYTICAL OR AZEOSTRINE HPL ASSTY, AFTER IT HAS BEEN FILTERED AGAIN THROUGH A 0.45 μm ACRODISC AND FILTER.

Recorded by: *John M. Kelly* Date: 04/28/06 Verified by: *Tom Bohman* Date: 10/25/07


To Page No. 119

AS PER K. KANE, THE REMAINING SAMPLE FROM 1044-117 WILL BE PASSED THROUGH A 0.45 μm FILTER AND SUBMITTED FOR ANALYTICAL TESTING.

AT THE REQUEST OF ANALYTICAL, AN ADDITIONAL SAMPLE OF FURTHERSON'S WITH SPERRY SUPERNATANT (FILTERED THROUGH 0.45 μm ACRODISC) WILL BE PROVIDED AND SUBMITTED FOR SIDE-BY-SIDE ANALYSIS/COMPARISON.

ONE ADDITIONAL SAMPLE FROM LOT # C184889 OF FURTHERSON'S REASONABLE NASH SPERRY, 50mg WAS OBTAINED. THE SAMPLE WAS STRIKEN TO SUSPEND THE MATERIAL, PLACED IN A SYRINGE @ ≈ 3300 RPM FOR 1 HOUR, THE SOLIDS DISINTEGRATED, AND THE SUPERNATANT POURED OFF.

THE SAMPLE WAS THEN PULSED INTO A 10cc SYRINGE AND FILTERED THROUGH SPIN, 2 μm , AND 1 μm ACRODISC FILTERS (SPIN). THE REMAINING SAMPLE WAS Poured INTO A NEW 10cc SYRINGE AND PLACED INTO A HOLE-IN-BACK CONSTANT PRESSURE SYSTEM AND ALLOWED TO SLOWLY FILTER THROUGH A 0.45 μm ACRODISC.

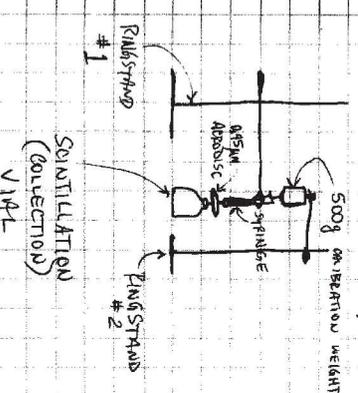

THE RESULTING SAMPLE WAS SUBMITTED TO ANALYTICAL (SPERRY) FOR TESTING, AND # 1044-119.

Recorded by: *John M. Kelly* Date: 04/28/06 Verified by: *Tom Bohman* Date: 10/25/07

To Page No. 119

Project No. _____

Book No. 1044

TITLE Soma SR 700 # 27-18-01c 04/28/06
PHYSICAL TESTING

From Page No. N/A

Physical Testing: Soma 700mg SR Tablets # 27-18-01c

Lot #	Target Wt. (mg)	Tablet #	Weight (g)	Thickness (mm)	Hardness (SCU)	Friability
27-18a-01c	1000	1	1.0144	6.57	28.6	Initial Wt (g) = 10.0540
		2	1.0024	6.56	26.7	Final Wt (g) = 10.0330
		3	1.0080	6.56	22.8	
		4	1.0087	6.53	27.8	
		5	1.0045	6.53	26.6	
27-18b-01c	1000	1	1.0065	6.68	28.3	Initial Wt (g) = 10.0200
		2	1.0063	6.70	25.6	Final Wt (g) = 9.9990
		3	1.0031	6.68	26.7	
		4	1.0008	6.70	28.3	
		5	0.9981	6.64	26.6	
27-18c-01c	1000	1	1.0049	6.80	28.2	Initial Wt (g) = 9.8440
		2	0.9922	6.73	26.2	Final Wt (g) = 9.8240
		3	0.9977	6.79	23.5	
		4	0.9815	6.76	27.2	
		5	0.9915	6.75	26.4	

Physical Testing: Soma 700mg SR Tablets # 27-19-01c

Lot #	Target Wt. (mg)	Tablet #	Weight (g)	Thickness (mm)	Hardness (SCU)	Friability
27-19a-01c	1000	1	0.9916	6.43	26.7	Initial Wt (g) = 9.9470
		2	0.9905	6.41	28.3	Final Wt (g) = 9.9220
		3	0.9813	6.43	23.4	
		4	0.9825	6.43	27.0	
		5	0.9861	6.48	25.1	
27-19b-01c	1000	1	1.0058	6.65	26.9	Initial Wt (g) = 10.0440
		2	1.0095	6.68	28.1	Final Wt (g) = 10.0260
		3	1.0051	6.63	27.2	
		4	1.0004	6.63	26.2	
		5	0.9996	6.62	27.2	
27-19c-01c	1000	1	0.9745	6.55	23.2	Initial Wt (g) = 9.8170
		2	0.9815	6.80	24.5	Final Wt (g) = 9.7990
		3	0.9845	6.58	23.1	
		4	0.9842	6.59	22.1	
		5	0.9714	6.55	20.8	

John's location of order job

ACTUAL TESTING PERFORMED BY: A. FAVARA.
ELECTRONIC COPY SUBMITTED TO B. JOHNS.
TABLETS ARE @ CLINICAL BATCHES, MANUFACTURED AT THE DECATUR FACILITY.

To Page No. N/A

Recorded by: [Signature] Date: 04/28/06 Verified by: QVA Date: 04/28/07

TITLE AZELASTINE SOLUBILITY EXPERIMENT - PEG 400

Project No. _____

Book No. 1044

05/01/06 121

From Page No. N/A

55.4500g = TARE WT. OF 100 FLINT GLASS JAR
 58.6771g = JAR + STIR ROD
 73.7000g = JAR + STIR ROD + ~15g PEG 400 (LOT # OP10-110)
 0.0075g = AZELASTINE @ 0.05% ADDED, ALLOWED TO MIX 15 MIN. STILL UNDISSOLVED. MIX ANOTHER 15 MIN., UNDISSOLVED. MIXED ANOTHER 30 MIN., SOLUTION IS CLEAR.

0.0076g = AZELASTINE ADDED TO BRING TO 0.1% AFTER 1 HOUR, SOLUTION IS CLEAR.

DECIDED TO RAISE AZELASTINE HCL TO 10 mg/mL, REFERENCING TABLE 5 OF THE NDA (pg. 03-63) WHICH STATES PEG 300 SOLUBILITY = 10.9g/mL

0.1500g = TARGET AZEL. HCL
 - 0.0151g = ALREADY ADDED.
 0.1349g = TO BE ADDED

0.1347g = AZELASTINE ADDED.

SAMPLE ALLOWED TO MIX OVERNIGHT.

DID NOT GO INTO SOLUTION, SAMPLE MILKY WHITE/OPAQUE.

~~N/A~~

To Page No. N/A

Recorded by: [Signature] Date: 05/01/06 Verified by: BN Babwani Date: 10/25/07

From Page No. 151

SPECTROPHOTOMETRY RESULTS:

12/05/06: FOLLOWING ABSENCE FROM WORK, SAMPLES WERE NOT TESTED IN AUGUST, AND ARE NOW PARTIALLY (AND SOME TOTALLY) EVAPORATED & WILL NOT BE TESTED.

(Signature)

~~NA~~

To Page No. 154

Recorded by:

John S. [Signature]

Date

12/05/06

Verified by:

Enl Bahwani

Date

01/25/07

From Page No. 154

- (A) 10:50a - ALL AZELASTINE HCL ADDED FIRST 0.1%
- 10:55 - ALL AZELASTINE HCL IN SOLUTION, SOLUTION CLEAR
- 11:00 - ALL NaCMC ADDED 0.2%
- 11:30 - SOLUTION CLEAR - ALL IS DISPERSED & IN SOLUTION, CLEAR

- (B) 10:00a - NaCMC ADDED 1.5%... CLEAR AFTER 25 MIN MIX
- 10:30 - ADDED PINCH OF AZELASTINE HCL, DOES NOT SEEM TO BE DISSOLVING
- 10:55 - ADDED ALL OF AZELASTINE HCL 0.1%, DOES NOT SEEM TO BE DISSOLVING
- 10:0p - AZELASTINE HCL FLOATING, NOT IN SOLUTION
- 3:00p - AZELASTINE HCL STILL UNDISPERSED / FLOATING IN SOLUTION AND SOME STUCK TO MIXING SHAFT.

EXPERIMENTS PERFORMED BY A. FAYARA.

~~NA~~

To Page No. 154

Recorded by:

John S. [Signature]

Date

12/05/06

Verified by:

Enl Bahwani

Date

01/25/07

(A.) 1025a - Add 40.15g Avicel RC-S91 (1.5%)
 mix 2:15 min.
 1045a - Add all AERASTRIVE HCL (0.1%)

RESULT IS A THICK, WHITE MIXTURE. AERASTRIVE HCL WEIGH
 GOES INTO SOLUTION. UNDISSOLVED CLUMPS ARE VISIBLE SUSPENDED
 IN A GEL-LIKE MIXTURE.

(B.) 1015a - AERASTRIVE HCL 0.1% added to water, mixed for
 5 minutes with dissolved,
 1050a - Add Avicel RC-S91 1.5% to the ^{aqueous} solution,
 mixed for 30 minutes. Mixture becomes a uniform
 solution.

SEPARATION occurs overnight. Top 1/3 is cloudy with some undissolved
 AERASTRIVE HCL observed. Bottom 2/3 is thick, white, opaque mixture
 that has settled to the bottom.

Sample sent to Amy Titch.

EXPERIMENT PERFORMED BY A. FARRERA.

~~NA~~

Recorded by: *[Signature]* Date: 12/05/06 Verified by: *[Signature]* Date: 10/25/07
 To Page No. 155

(A.) AERASTRIVE HCL 0.15% added to water, mixed for 5 minutes with
 dissolved and clear.
 Avicel RC-S91, 1.5% added next, mixed for 30 minutes, does
 not mix evenly.
 No gel formation is observed.
 While the sample settles to bottom and is opaque white.
 Top half is cloudy, but TRANSLUCENT.

(B.) Avicel RC-S91, 1.5%, added to water
 mixed for 15 minutes, resulting in a thick, opaque mixture.
 A sample of this is submitted to Amy Titch.

EXPERIMENT PERFORMED BY A. FARRERA.

~~NA~~

Recorded by: *[Signature]* Date: 12/05/06 Verified by: *[Signature]* Date: 10/25/07
 To Page No. 154

Project No. _____

Book No. 1044

TITLE METHOCEL SOLUTION EXPERIMENTS

02/12/07

From Page No. 1048-184... EXPERIMENTS BEGUN BY A. FAUVER

Viscosity Measurement of Methocel Solutions

① METHOCEL ES (Lot # TB0901240) @ 5%

MEASURE 50 mL P. WATER, HEAT TO 80°C. ADD METHOCEL (5g) AND STIR UNTIL DISPERSED. ADD 50 mL OF COLD P. WATER. CONTINUE TO MIX UNTIL SOLUTION IS CLEAR. ALLOW TO REACH ROOM TEMPERATURE. MEASURE VISCOSITY USING BROOK FIELDS VISCOSITOMETER

Visc. Reading: $0.2 \times \text{FACTOR} 100 = 20 \text{ cps}$

② METHOCEL ES (Lot # TB0901240) @ 25%

CREATED WITH SIMILAR PROCEDURE, EXCEPT QUANTITIES DOUBLED FOR RATE OF ADDITION.

Visc. Reading: $3.0 \times \text{FACTOR} 100 = 300 \text{ cps}$

③ METHOCEL K100LVCR (Lot # TJ2012N31) @ 2.5%

CREATED WITH SIMILAR PROCEDURE TO #2
Visc. Reading: $7.6 \times 100 = 760 \text{ cps}$

④ METHOCEL K100LVCR (Lot # TJ2012N31) @ 10%

CREATED WITH SAME PROCEDURE
Viscosity Reading: ERROR :: VISC > 10,000 cps

⑤ METHOCEL K100MCR (Lot # UD11012N32) @ 2.5%

CREATED WITH SAME PROCEDURE
Viscosity Reading: ERROR :: VISC > 10,000 cps

⑥ METHOCEL K100MCR (Lot # UD11012N32) @ 10%

COULD NOT BE PRODUCED.

To Page No. _____

Recorded by: *[Signature]*

Date

02/13/07

Verified by:

[Signature]

Date

24
02/07

TITLE

ASTELIN 03-33 PUMP CLEANING EVAL.

Project No. _____

Book No. 1044

03/01/07

157

From Page No. N/A

REPORT No: PPD-07-07P VERSION: 00

PURPOSE: TO DETERMINE IF THE CURRENT PUMP CLEANING PROCEDURE PROVIDED TO THE CONSUMER WITH COMMERCIAL ASTELIN IS SUITABLE FOR FORMULATION # 03-33.

PACKAGING: BATCH 03-33-03c, 30 mL-FILL HDPE V-BOTTOM WITH VALOIS SPLIT RING PUMP. (34.5 mL, HDPE, CODE # G1834, LOT # 000002836).

EXPERIMENT: REMOVE DUST CAPS & CAPS FROM 20 BOTTLES OF AZELASTINE HCL SOLUTION 0.1% w/v, FORM. 03-33, BATCH 03-33-03c. IN A GOOD ACTUATE EACH PUMP AT LEAST 6 TIMES OR UNTIL A FINE SPRAY APPEARS. LABEL HALF THE BOTTLES "UPRIGHT" AND THE OTHER HALF "PRONE." NUMBER EACH SET 1 THROUGH 10. STORE THE TEN BOTTLES UPRIGHT IN A 40°C/AMBIENT RH CHAMBER. STORE THE TEN BOTTLES PRONE IN THE SAME CHAMBER. AFTER 2 WEEKS STORAGE, REMOVE BOTTLES AND ATTEMPT TO ACTUATE EACH.

OBSERVATIONS: ALL 20 SAMPLES PERFORMED IN SIMILAR FASHION. NOTHING NOTABLE OBSERVED DURING ^{THE} SPRAYING. ALL SAMPLES PERFORMED NORMALLY. INITIAL ^{INITIAL} *[Signature]*

SAMPLES WERE RE-CAPPED, NECK-CLIPS WERE REPLACED, AND ALL SAMPLES WERE PLACED IN 40°C/AMBIENT STORAGE CHAMBER AQ IN ROOM #142.

THESE SAMPLES WILL BE RETRIEVED AFTER 2 WEEKS AND EACH WILL BE ACTUATED (OR AT LEAST ATTEMPTED). IF THE PUMPS WILL NOT ACTUATE, THE PATIENT INSTRUCTIONS FROM THE COMMERCIAL PRODUCT WILL BE FOLLOWED.

To Page No. 158

Recorded by: *[Signature]*

Date

03/05/07

Verified by:

[Signature]

Date

3/26/07