

www.archive.org 415.561.6767 415.840-0391 e-fax

Internet Archive 300 Funston Avenue San Francisco, CA 94118

#### AFFIDAVIT OF ELIZABETH ROSENBERG

- 1. I am a Records Request Processor at the Internet Archive, located in San Francisco, California. I make this declaration of my own personal knowledge.
- 2. The Internet Archive is a website that provides access to a digital library of Internet sites and other cultural artifacts in digital form. Like a paper library, we provide free access to researchers, historians, scholars, and the general public. The Internet Archive has partnered with and receives support from various institutions, including the Library of Congress.
- 3. The Internet Archive has created a service known as the Wayback Machine. The Wayback Machine makes it possible to browse more than 450 billion pages stored in the Internet Archive's web archive. Visitors to the Wayback Machine can search archives by URL (i.e., a website address). If archived records for a URL are available, the visitor will be presented with a display of available dates. The visitor may select one of those dates, and begin browsing an archived version of the Web. Links on archived files in the Wayback Machine point to other archived files (whether HTML pages or other file types), if any are found for the URL indicated by a given link. For instance, the Wayback Machine is designed such that when a visitor clicks on a hyperlink on an archived page that points to another URL, the visitor will be served the archived file found for the hyperlink's URL with the closest available date to the initial file containing the hyperlink.
- 4. The archived data made viewable and browseable by the Wayback Machine is obtained by use of web archiving software that automatically stores copies of files available via the Internet, each file preserved as it existed at a particular point in time.
- 5. The Internet Archive assigns a URL on its site to the archived files in the format http://web.archive.org/web/[Year in yyyy][Month in mm][Day in dd][Time code in hh:mm:ss]/[Archived URL] aka an "extended URL". Thus, the extended URL http://web.archive.org/web/19970126045828/http://www.archive.org/ would be the URL for the record of the Internet Archive home page HTML file (http://www.archive.org/) archived on January 26, 1997 at 4:58 a.m. and 28 seconds (1997/01/26 at 04:58:28). A web browser may be set such that a printout from it will display the URL of a web page in the printout's footer. The date indicated by an extended URL applies to a preserved instance of a file for a given URL, but not necessarily to any other files linked therein. Thus, in the case of a page constituted by a primary HTML file and other separate files (e.g., files with images, audio, multimedia, design elements, or other embedded content) linked within that primary HTML file, the primary HTML file and the other files will each have their own respective extended URLs and may not have been archived on the same dates.
- 6. Attached hereto as Exhibit A are true and accurate copies of printouts of screenshots of the Internet Archive's records of the archived files for the URLs and the dates specified in the attached coversheet of each printout.



7. I declare under penalty of perjury that the foregoing is true and correct.

DATE: 11/02/2020


Elizabeth Rosenberg

## **JURAT**

State/Commonwealthof	TEXAS )			
	)			
LCity ✓County of	Wise )			
On 11/02/2020 , bet	fore me, Shannara Ann Franzel			
Date	Notary Name			
the foregoing instrument w	as subscribed and sworn (or affirmed) before me by:			
	Elizabeth December			
	Elizabeth Rosenberg  Name of Affiant(s)			
	Name of Amanus)			
<ul><li>Personally known to me - 0</li></ul>	DR			
D. Duessad to use on the besie of	falso and sof			
Proved to me on the basis of the oath of				
	f satisfactory evidence: passport			
	Type of ID Presented			
minimining.	7			
Shannara Ann Franzel				
ID NUMBER	WITNESS my hand and official seal.			
132662714 COMMISSION EXPIRES				
August 27, 2024	Notary Public Signature: Nannara Alan Franzel			
	Channeys Ann Evental			
Notary Name: Shannara Ann Franzel				
Notary Commission Number: 132662714				
Notary Commission Expires: 08/27/2024				
	Notarized online using audio-video communication			
DESCRIPTION OF ATTACHED	DOCUMENT			
Title or Type of Document: Aff	idavit of Authenticity			
Document Date:11/02/2	2020			
Number of Pages (including not	tarial certificate):3			

# Exhibit A

nttps://web.archive.org/web/20170705124424/https://ryortho.com/2016/07/th e-eleven-best-sports-medicine-technologies-for-2016/	

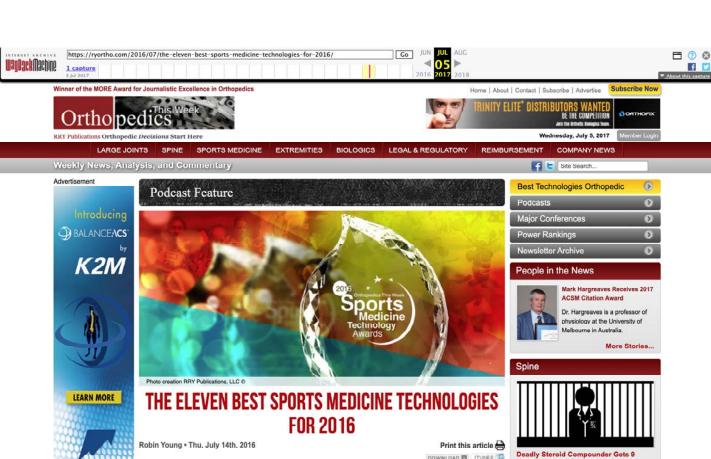

## REDEFINING





Come visit us NASS Summer Meeting

ORTHOFIX

Advertisement

Mobi-C® **Cervical Disc** 



Learn more about Mobi-C

DOWNLOAD TONES The 11 winning companies with the BEST new sports medicine technologies for 2016 are: AlloSource, Arch Day Design (Smith & Nephew). CyMedica Orthopedics. Lipogems. Micro-Imaging Solutions. Orthopedic Sciences.

#### Rewarding Innovation and Perspiration

Össur, Stryker, Topical Gear, VICIS, Vomaris Innovations and ZipLine Medical.

This annual 2016 award rewards inventors, engineering teams, physicians and their companies who've created the most innovative, enduring and practical products and technologies to treat injured athletes in ANY activity—whether sports, industry, military or healthcare. To win the Orthopedics This Week Best New Technology Award for Sports Medicine care, a new technology must meet the following criteria:

- 2. Have long term significance to the problem of treating the injured athlete. Does this technology have staying
- 3. Solve a clinical problem. To what extent does this technology solve a current clinical problem or problem that is inadequately solved today?
- 4. Does it have the potential to improve standard of care?
- 5. Is it cost effective?
- 6. I would use it.

Our panel of physicians score every submission on a scale of 1 to 5 (5 being the highest score) for each of the above

We and our panel of surgeons were impressed that inventors—despite ever growing hurdles to innovation and entrepreneurism in medicine—still managed to create a solid group of nearly 40 new products to submit for the 2016 Orthopedics This Week Sports Medicine Technology Awards.

We offer our thanks and deep appreciation to the engineering teams, surgeon inventors and the following companies for submitting their best ideas this year:

ConMed

Aperion Biologics, Inc. Active Implants LLC AlloSource® Augmenta Spine Arch Day Design / Smith & Nephew Arthrex, Inc.

CollPlant LTD



#### The Judges

Our intrepid panel of surgeon judges were:

· James Andrews, M.D.: Dr. Andrews is the Dean of Sports Medicine and Clinical Professor of Orthopaedic









Patent Owner Cellect, LLC Ex. 2102, p. 6 Leading the way for both one and two-level cervical disc arthroplasty

ZIMMER BIOMET

SI-BONE

iFuse-3D™

A New Dimension in SI Joint Fusion<sup>®</sup>

> Click Here to learn more

Surface Matters. xperience Matters





#### All the Benefits of Ti-Bo

- Over 10,000 Devices Implanted · 3 Years' Clinical Experience
- Roughened Surface
- us Titanium
- Now with a Warranty\* Aga
- Pseudarthrosis

SPINAL Click Through to Learn More...

\* Now Warrantied **以下是一个的** 

Surgery at the University of Alabama Birmingham Medical School, the University of Virginia School of Medicine, the University of Kentucky Medical Center, and the University of South Carolina Medical School. Dr. Andrews also holds a Doctor of Laws Degree from Livingston University, Doctor of Science Degree from Troy State University and a Doctor of Science Degree from Louisiana State University. He is the founding partner and medical director of the institute that bears his name—The Andrews Institute in Guif Breeze, Florida. Dr. Andrews is mentor for more than 150 orthopedic/sports medicine fellows and more than 30 primary care sports medicine fellows. He has made major presentations on every continent and has authored numerous scientific articles and books. Dr. Andrews is a ard of Orthopaedic Surgery and the American Academy of Orthopaedic Surgeons. He has served on the Board of Directors of the American Orthopaedic Society of Sports Medicine and the Arthroscopy Association of North America and the International Knee Society. At present, Dr. Andrews serves as Co-Medical Director for Intercollegiate Sports at Auburn University, Senior Orthopaedic Consultant for Intercollegiate Athletics at the University of Alabama and consultant for the athletic teams of Troy University, University of West Alabama, Tuskegee University and Grambling University.

- Brian Cole, M.D.: Cole is Associate Chairman and Professor, Department of Orthopedics, Chairman, of the Department of Surgery at Rush and Section Head, Cartilage Restoration Center at Rush. Dr. Cole is also a Professor in the Department of Orthopedics with a conjoint appointment in the Department of Anatomy and Cell Biology at Rush University Medical Center in Chicago, Illinois. In 2011, he was appointed as Chairman of Surgery at Rush Oak Park Hospital. He also serves as the head of the Orthopedic Master's Program and trains resident: and fellows in sports medicine and research. He lectures nationally and internationally, and through his basic science and clinical research has developed several innovative techniques for the treatment of shoulder, elbow and knee conditions. He has published more than 1, 000 articles and has published 5 widely read textbooks in orthopedics. In 2009, Dr. Cole was selected as the "NBA Team Physician of the Year." He is the team physician for the Chicago Bulls NBA Basketball team, co-team physician for the Chicago White Sox Major League Baseball team and DePaul University in Chicago
- John Dorris, M.D.: Dr. Dorris is one of busiest sports and general orthopedic physicians in the United States. A member of the American Academy of Orthopedic Surgeons, Clinical Orthopedic Society, American Association of Sports Medicine, Southern Orthopedic Association, and Georgia Orthopedic Society, Dr. Dorris has built a strong practice treating joint pain, sports injuries (ACL tear, meniscus tear, rotator cuff tear, tennis elbow, golfer's elbow, labral tear), tendon pain (elbow pain, heel pain, knee pain), fractures, and other orthopedic issues. His medical education is from Tulane University, internship and residence at Emory University and fellowship at Steadman Hawkins clinic in Vail, Colorado.
- David Petron, M.D.: Petron is Assistant Professor (clinical) at the University of Utah and is a primary care orthopedic/sports medicine specialist. He originally trained in family practice and then completed a fellowship in primary care orthopedics and sports medicine at Michigan State University. Dr. Petron is the Director of Primary Care Sports Medicine and the team physician for the University of Utah athletic teams. He is also a team physician for the U.S. Ski Team and is the Medical Director for the Park City Ski Injury Clinic. Dr. Petro specializes in non-operative treatment of the spine, shoulder, knee, ankle, and foot problems. He also has a special interest in overuse injuries that lead to tendonitis. He utilizes the latest treatments for this, including shock wave therapy (Sonarex), for chronic tendon problems. Dr. Petron also has a special interest in all non-operative ents for osteoarthritis and performs various injections (such as Synvisc) to alleviate joint pain
- Tom Vangsness, M.D.: Dr. Vangsness is a Professor in the Department of Orthopaedic Surgery at the Keck School of Medicine of University of Southern California (USC). Dr. Vangsness is internationally recognized for his expertise in shoulder and knee surgery and has conducted and published extensive research with a focus on such topics as shoulder and ligament biomechanics, molecular biology and sports medicine, allograft transplantation, radiofrequency energy effects on articular cartilage, meniscal cartilage and ligaments. Dr. Vangsness is a team physician for USC Athletics and has a special interest in sports injuries. A recipient of numerous university and national academic awards as a teacher and researcher, Dr. Vangsness is also on the editorial board of the .lourn. of Arthroscopy and The Journal of Knee Surgery. He is a member of several other prominent organizations including the American Orthopaedic Society for Sports Medicine, American Shoulder and Elbow Surgeons, and The International Society of Orthopaedic Surgery and Traumatology. Dr. Vangsness has recently published a book entitled "The New Science of Overcoming Arthritis

Here, then, are the 11 best new sports medicine technologies for 2016 arranged in alphabetical order.

#### AlloSource®

Winning Technology: ProChondrix® Cartilage Restoration

Inventors: Carolyn Barrett, MBA - Scientist

Engineers: Matt Southard - R&D Engineer, Meghan Gabriel -

Manufacturing Engineer

(L to R): Robin Young, Peter Stevens, Ph.D., MBA, Susan Tice, ProChondrix is a cellular three dimensional fresh allograft cartilage MBA, Carolyn Barrett, BS, MBA, Meghan Gabriel, Tom Cycyota natrix. This matrix provides live, functional cells and other

biological components necessary for the natural repair and regeneration of damaged cartilage tissues.

The implant also has native growth factors, an extracellular matrix which helps support a structure for cellular migration and adhesion and viable chondrocytes. It is a purely natural solution intended to maintain its homologous use from donor to recipient.

ProChondrix naturally contains the three key components (scaffold, signal, and cells) necessary for the repair and regeneration of damaged cartilage tissues. It naturally matches biomechanical and biochemical properties of norma hyaline cartilage, may restore a smooth articular cartilage surface, relieve patient symptoms and improve function.

ProChondrix is a cost-effective, single-stage alternative to marrow stimulation procedures, chondroplasty or fullthickness osteochondral replacements. It is implanted arthroscopically or in a minimally invasive procedure. Finally, ProChondrix eliminates the need for autografting, which risks donor-site morbidity.

#### Arch Day Design and Smith & Nephew

Winning Technology: AccuPass Direct Suture Passer

Inventors: Roger Pisarnwongs, Tom Weisel Seth Foerster

Engineering Team: Roger Pisarnwongs, Tom Weisel, Seth Foerster, Vince Tangherlini, Greg Bain

The best ideas are the simplest ideas. One of the challenges of passing suture through tissue during an arthroscopic procedure is Woelfel, Danny Moriarty, David Miller, Tyler Linna the difficulty to grab a small loop of suture just barely peeking out

Smith & Nephew (L to R): Gave Veloz, Eric Duggan, Stepi

or lost among 'fluffy' tissue. The AccuPass Suture Passer presents the surgeon with a large loop of suture after passing so the user can easily grab the suture and pull it out the cannula to complete the pass. This creative mechanism is flexible enough to allow usage within a variety of standard needle shapes for extensive anatomical access.

Click Here

www.si-bone.com/risks for safety information

#### Sports Medicine



Which is Safer, Organized Sports or Free





#### Large Joints



**Exactech Wins Japanese Approval** for Joint Replacement Devices



Free Surgeries Honor Mandela Day







Patent Owner Cellect, LLC Ex. 2102, p. 7 Inis solution works well where joint spaces are very tight and don't allow for aggressively free needle movement. With the AccuPass you can deploy a loop of suture on the far end of the tissue so large that when the surgeon removes the needle it does not appreciably pull back and affect the loop size. That large loop is easy for the user to grab and pull out the cannula with this simple, intuitive, one-handed device. AccuPass uses a specially designed Nitinol wire that hooks the suture but also pulls it deep into the needle tube. Once the needle pierces the tissue, the unique wire design ejects the captured suture out the tube to become a large loop on the far side. Even if the loop is pulled back slightly during removal of the needle there is still a large target (suture loop) remaining that is easy for the surgeon to grab and retrieve. The wire form has been



Arch Day Designs (L to R). Tom Weisel and Roger Pisarnwongs

optimized to handle a variety of suture sizes and to be used several times in a case. It has also been designed for inexpensive manufacturing that allows it to be built in to a price competitive single patient use device.

#### CyMedica Orthopedics, Inc.

#### Winning Technology: e-Vive with CyMotion Technology

Inventore: Struan Coleman, M.D., Ph.D., David Saar, Calvin Domenico, Lee Knox, Edison Gieswein, Marlina Kessler, Joshua Butters

Engineering Team: Joshua Butters, Marlina Kessler, David Saar, Benjamin Ciguere, Lee Knox, Dharmavirsinh Jhala, Shane McCauley



(L to R): Bryti Ketchum, Rob Morocco, Struan Coleman, M.D., Ph.D., Robin Young, Ron Yagoda, Chairman, Josh Butters

CyMedica's e-Vire ™ muscle activation system with CyMotion™ neuromuscular electrical stimulation (NMES) technology delivers an individualized, home based, wireless, app-controlled NMES therapy. The system also includes a telehealth data capture and communications capability to transfer real time data to physicians through the cloud based intelliHUB™ provider portal.

The primary objectives of sports related knee ligament surgery is to restore knee function to pre-injury levels and promote long-term joint health. Quadriceps strength and endurance are of vital importance for normal knee joint function, so restoring normal quadriceps function after knee joint injuries is an essential component of rehabilitation. Typically rehab programs to strengthen quadriceps focus on voluntary exercises which have hit or miss outcomes. Home Neuromuscular electrical stimulation (NMES) systems are a promising intervention to use after ligament repair removing the barrier of muscle inhibition on behalf of the patient and activating the quadriceps muscle.

The e-Vive system when combined with CyMotion changes the standard of care for ACL quadriceps strengthening through the remote patient monitoring capabilities and the precise control features of the wirelessly sensored e-Vive app and garment. Clinicians can engage patients in their rehabilitation protocols through interactive goal setting of total range of motion, pain control measures and tracking total steps while also receiving data via the real time telerehab evaluation.

#### Lipogems

#### Winning Technology: Lipogems System

Inventors: Dr. Carlo Tremolada

The Lipogems system is a closed-looped point-of-care device that gently washes and mechanically micro-fragments autologous adipose tissue. This tissue can be used for cushioning, repairing, and restructuring soft tissue defects. For example, under ultrasound guidance, the Lipogems processed tissue can be injected directly into the damaged native tissue to provide a structural graft that sticks to the area and provides a regenerative scaffold.

(L to R): Mollie Pett, Carl Llewellyn, Robin Young

Lipogems is an FDA compliant option for utilizing the well-known regenerative properties of autologous adipose tissue. Because Lipogems maintain the structural characteristics of the tissue, they are not considered more than minimally manipulated and fall under the CFR1271.361 exemptions, enabling physicians to use the device point-of-care under current FDA regulations.

The viscous properties of the tissue graft allow Lipogems to be injected through a small gauge needle, and the intact stromal vascular niche provides a naturally derived regenerative matrix. The adhesive quality of adipocytes allows the Lipogems to remain where they are injected, providing structural support while remodeling of the tissue occurs. Adipose tissue contains cells and tissue that are capable of regeneration and are optimized by the washing and resizing of the Lipogems. This tissue may exhibit anti-microbial and immune-modulatory characteristics and natural anti-inflammatory effects that may provide a long lasting and natural alternative to steroids, and potentially a regenerative alternative or adjunct to surgical procedures.

#### Micro-Imaging Solutions, LLC

Winning Technology: Distally mounted CMOS Sensor for disposable/reusable and rigid/flexible scopes

Inventors: Jeffrey Adair

This technology allows for extremely small and higher resolution CMOS Sensors to be placed in distal tip of endoscopes. Current cameras and endoscopes (all types) have remained unchanged for decades and degrade from "Day 1" of utilization thus causing undesirably compromise of visual clarity at an extremely inconvenient time! This technology allows for flexible or rigid





Advertisement



Advertisement



Learn more at osseus.com

#### Most Popular

\*\*\*\*

Trousdale v. MacDonald: The Cemented All Poly Tibia in the Active <60 Patient



Alejandro Badia, M.D. Named Finalist in Health Care Heroes Awards



Sports Medicine's Role in the Opioid Epidemic



Stryker Crosses 100,000 Robotic Procedure Milestone



How to Cut THA Stay From 3 to 0



Patent Owner Cellect, LLC Ex. 2102, p. 8

applications as a reusable or a disposable with significant value. Utilization of this technology can eliminate the five things that degrade current and problematic technology, it has the real benefit (L to R): Robin Young, Eric Dahlinger of reducing repair costs as well as OR "down time" and the need to

acquire costly surplus inventory to avoid catastrophes. Every patient can get a new and sterile scope and new applications for scopes are achievable-i.e., flexible HIP Arthroscopes. It clearly can reduce the true cost of ownership and provide OR efficiencies not currently available.

#### Orthopedic Sciences, Inc.

Winning Technology: Piranha Arthroscopic Shaver, Talon Burr and Harrier

Inventor: James K. Brannon, M.D., FAAOS

The Piranha is an arthroscopic shaver with improved fluid flow in single axis rotation. Arthroscopic shavers traditionally oscillate. Many orthopedists would prefer a single axis (direction) of rotation for a smoother surface but, unfortunately, fluid flow with the single axis rotation is diminished. Thus, oscillation is the most frequently used handpiece mode. The Piranha is a standard arthroscopic shaver with improved fluid flow characteristics. The Piranha is available in full radius and fast cut configuration. The lengths are 132 mm and 180 mm. Diameters are 4.5 short and long, and 132 mm for a 3.5 mm



James K. Brannon, MD, FAAOS

Piranha delivers PRECISE, improved and level tissue debridement that is thorough. Many basic tools for arthroscopy are outdated and more thorough debridement can be achieved when visibility is enhanced, shaver fluid-flow is improved, and burrs are more effective at contouring articular surfaces.

The Talon Burr and the Harrier Burr are disposable, single use, unshrouded handheld arthroscopic burrs. Both fit onto Smith & Nephew, Stryker, and Linvatec handpieces. The Talon Burr is particularly useful at the patella femoral sulcus region. The Harrier Burr is tailored for more surface congruent chondroplasty, i.e., a microchondroplasty and improved cartilage contouring during arthroscopy of the hip, shoulder, and knee joints. The Harrier Burr rapidly contours the cartilage surface more effectively and rapidly than an oscillating arthroscopic shaver or a standard burr. Therefore, more congruent articular surfaces can be realized postoperatively in the performance of a microchondroplasty.

#### Össur

#### Winning Technology: Rebound® PCL

Inventors: Arni Thor Ingimundarson, Harry Duane Romo, Bjorn Omarsson, Sri Chetlapali, Thorvaldur Ingvarsson, M.D., Ph.D. and Robert Francis LaPrade, M.D., Ph.D.

Engineer: Jared Olivo, Jane Lee, Bjorn Omarssom, Sri Chetlapali

The Rebound PCL is the world's first dynamic force PCL brace The brace is designed to apply a physiologically correct dynamic force, optimum for rehabilitation of posterior cruciate ligar (PCL) ruptures, whether during functional (non-surgical) treatment (I. to R): Bret Carter, Robin Young or post-surgical reconstruction. The dynamic load of the Rebound



PCL on the tibia is generated by applying an anterior directed dynamic force on the calf area, and opposing counter forces on the anterior aspect of the leg.

The Rebound PCL is indicated for rehabilitation of posterior cruciate ligament (PCL) ruptures, whether post-surgical reconstruction or during functional (non-surgical) treatment.

In a PCL deficient knee, gravity and the forces on the joint from the hamstring muscles can potentially cause the tibia to be positioned in a posterior subluxed location relative to the femur (posterior sag). Healing of the PCL in an elongated position can lead to chronic instability and disability. Therefore, anatomically correct positioning during conservative or surgical rehabilitation is important.

Current PCL braces either act as a barrier to normal motion or apply a static load. The Rebound PCL brace, by contrast, applies physiologically correct dynamic force to the PCL during rehabilitation. The Dynamic Tension System™ (DTS) increases the load on the tibia as the knee goes through flexion, providing an Anterior Drawer force that reduces the load on the PCL, allowing it to heal. The Rebound PCL allows for specific load adaptation according to the patient's individual anatomy and rehabilitation related requirements. For complex ligament reconstruction with PCL injury, the Rebound PCL applies physiologically correct dynamic support for PCL and ensures medial/lateral

#### Stryker Corporation


### Winning Technology: SonicAnchor

R&D Team: Dr. Bernd Robioneck, Dr. Robin Buescher, Joerg Thon Helge Giersch, Rouven Borowsky, Salvador Jimenez, André Wilke

Marketing Team: Torn Popeck, Boon Quah, Mark Baldisserotto,

SonicAnchor is a truly disruptive new suture anchor design. It uses (L to R): Zachary Meli, Liam McCarthy, Robin Young ultrasound to liquefy a polymeric material which interdigitates with surrounding cancellous bone and securely fixes suture and the soft tissue

The SonicAnchor system is a small size implant made of a PLDLLA co-polymer that is fully bioresorbable. After a suture is connected, ultrasound is applied, the implant liquefies and interdigitates with the surrounding cancellous

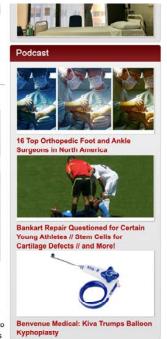

# Alignment NUVASIVE



Biomechanical tests have shown that SonicAnchor provides superior fixation when compared to competitive suture anchors of similar size. With a diameter of 2.5mm, the SonicAnchor needs little surface real estate and no cortical bone for superior for fixation. The implant is radiolucent and MRI safe

The SonicAnchor was developed to provide superior pull-out and fatigue strength even in de-corticized bone. Median fatigue limit of the SonicAnchor vs. a conventional suture anchor is approximately +95%. The static pull difference is even higher

SonicAnchor uses PLDLLA which has not been shown to create any adverse events in a four-year clinical trial.

Compared to a threaded interface, the liquefied material results in a larger interface area and thereby provides a higher bonding. As a consequence, implants with a smaller footprint can be developed. In addition, the technology allows using next generation bioresorbable polymers so that hardware removal of the implant is not required and MRI safety is given.

#### **Topical Gear**

Winning Technology: ACL Tube

Inventors: Bill Bue

Engineers: Betsy Russel, Bob McCune

The ACL Tube is a form-fitting thigh sleeve with two embedded buttresses that applies focused topical pressure over the medial quadriceps and hamstrings muscles. Topical pressure stimulates (L to R): Robin Young, Bill Blu, Michael Decker, Ph.D. sensory receptors located in these muscles and enhances the



central nervous system's ability to detect, interpret and utilize sensory information to produce a coordinated muscular response for dynamic knee joint stability. This neuromuscular response is learned and retained when wearing the ACL Tubes during training or rehabilitation.

Proprioceptive stimulation to the muscles elicits a high fidelity flow of sensory information for the CNS to detect, interpret and formulate quality movement ultimately circumventing neuromuscular deficits. This technology can be used to augment training and rehabilitation for the prevention or functional recovery of musculoskeletal issues. More importantly, it can be used in combination with other devices (i.e., knee bracing) to optimize their positive effects (i.e., protection) while minimizing their negative effects (i.e., atrophy). Hence, the ACL Tube provides a robust clinical tool that improves upon current standards of care.

#### VICIS, Inc

#### Winning Technology: ZERO1

Inventors: Dave Marver, Samuel Browd, M.D., Ph.D., Per Reinhall Dh D

Engineering Team: Samuel Browd, M.D., Ph.D., Per Reinhall,

VICIS has developed a completely redesigned football helmet named the ZERO1 which mitigates linear and rotational impact forces to the head. The helmet utilizes four discreet layers of materials to reduce injuries as the result of impacts to the head:



(L to R): Per Reinhall, Ph.D., Dave Marver, Samuel Browd, M.D., Ph.D

- 1. LODE Shell: The outer layer that absorbs impact load by locally deforming, like a car bumper.
- 2. RFLX Layer: A columnar structure that moves omni-directionally to reduce linear and rotational forces. The LODE Shell and RFLX Layer work together to reduce impact forces.
- 3. ARCH Shell: VICIS' AXIS Fit System determines the relationship between head length and breadth measurements (commonly called an "aspect ratio") for each player.
- 4. FORM Liner: The FORM Liner is conforms to the player's unique head topography and distributes pressure uniformly around the head

Concussions are a growing issue in all levels of sport, with the most attention currently being paid to football at all levels of the game.

VICIS' ZERO1 uses a flexible outer shell over a layer of flexible columns, so when a blow is delivered, the outer layers deform, absorbing energy, then spring back while a stiff inner shell cradies the player's head. Because the flexible columns can move in all directions, the system is better at buffering off-center blows that cause the brain to rotate inside the skull, mitigating acceleration while still preventing catastrophic head trauma.

The ZERO1 will be commercially available by the second half of 2016, and has received commitments to use the helmet from several teams at every level of football.

#### Vomaris Innovations, Inc.

Winning Technology: Procellera® Antimicrobial Wound Dressing with Advanced Microcurrent Technology®

Inventors: Jeffry Skiba, Larry Schneider

Engineers: Troy Paluszcyk, Joseph Del Rossi



(L to R): Troy Palusczyk, Michel Tan

Advanced Microcurrent Technology® (AMT) is a dressing for wound and inclaims a contract of the and incisional care which employs embedded microcell batteries

which generate microcurrents in the presence of moisture. Procellera's goal is to mimic the body's own physiologic electrical signals to direct cell migration and re-epithelialization. It has demonstrated efficacy against a br spectrum of microbes including drug-resistant and biofilm-forming pathogens. This unique wound dressing creates an optimal environment for wound healing.

AMT is distributed as JumpStart<sup>TM</sup> Antimicrobial Wound Dressings exclusively by Arthrex, Inc. in the United States.

Over the past couple years there have been some very public reports on locker room-, playing field or equipmentacquired infections with antibiotic resistant microbes that pose significant threats to athletes.

AMT is a *non-pharmacologic* approach to broad-spectrum antimicrobial efficacy and is available in configurations designed specifically for the unique needs of orthopedic surgeons, team doctors/trainers, and their patients. AMT also stays in place and moves with joints, allowing patient/athletes unrestricted range of motion (ROM) and mobility so they can return to activity.

#### **ZipLine Medical**

#### Winning Technology: Zip Surgical Skin Closure

Inventors: Amir Belson, M.D., Eric Storne, MBA, Eric T. Johnson, Robert R. Ragland, Phillip C. Burke, Luke Clauson

Engineering Team: Kei Ichiryu, MS, Zach Kimura, Alan Schaer. Daren Stewart, MS, Melissa Guerrero, Jeremy Edinger, Emily Cullinan, Lori Munoz, Julie Ridgeway, John Tighe, Eric Storne



The Zip® Surgical Skin Closure provides a non-invasive alternative (L to R): Robin Young, Eric Dahlinger to staples, sutures and glue for surgery and lacerations. It is used

by surgeons (orthopedic, general, cardiac, ob/gyn, plastic/reconstructive, dermatology) and emergency department physicians, and because of the device's ease of use, no suturing skill is required so the closure task using the Zip device can be delegated to a Physician's Assistant or RN.

This Zip is applied to intact skin immediately surrounding a surgical incision, and is placed at the end of the procedure, when the skin layer is to be closed. It replaces conventional means of skin closure such as staples, sutures and glue. It utilizes a strong yet skin-friendly hydrocolloid adhesive that lasts for 14 days, and is easily removed. The incision is closed by sequentially tightening a series of ratcheting straps, similar to "zip-ties," to produce the decired incision closure and tension.

The Zip consists of biocompatible polymeric materials, including hydrocolloid adhesive, polyurethane monofilm, nylon injection-molded tensioning mechanism and polyester force-distributing plates. It is a single-use device and is provided sterile in a sealed pouch suitable for aseptic transfer to a sterile field.

Leave a Reply		
Name		
	CAST NOTE OF THE STATE OF THE S	
Email Address (will not be pul	olished)	
Website		
Comment:		

LARGE JOINTS SPINE SPORTS MEDICINE EXTREMITIES BIOLOGICS LEGAL & REGULATORY REIMBURSEMENT COMPANY NEWS



Copyright 2017 Orthopedics This Week Inc. All rights reserved. Sitemap | Privacy Policy Become a Member | Advertise | Login | Contact Website design by Inverse Paradox